

200203 - PdAg + 550 K O2 + H2

Monday, February 3, 2020 9:39 AM

Plan

Test reproducibility of 550 K 5 min oxidation and establish protocol

Logbook

Legend - Plan in purple

Flash anneal to 800 K

Set up water cooling

Open gate valve

Measure clean Ag and turn off O3 (8:50 AM)

XPS: clean Ag - 66 degrees

Switch on evap

Connect cooling

9:45 AM

- **Pd evaporation:** 2.50 A, 1000 V, 10.0 mA, 1.7-2.0 nA for **1 min 20 sec** at 300K

Goal: 0.10 ML

$0.14/1.5 = 0.0933$

$0.1/0.0933 = 1.0718$

$1.0718 * 2.0/1.7 = 1.2609$

$0.26 * 60 = 15.6$

Switch on X-Ray source before transfer

XPS: Pd deposition - 0 degrees

- Pd area: 11 kcps eV
- Ag area: 713 kcps eV
- Pd/Ag: $11/713 = 0.0154$ --> **0.09(3) ML GREAT**

XPS: Pd deposition - 66 degrees


- Pd area: 11 kcps eV
- Ag area: 275 kcps eV
- Pd/Ag: $11/275 = 0.04$ --> **0.10(8) ML TOP**

Done at 10:15 AM

Turn up O2 at 10:30 AM

- **Oxidized** in 1.92(7) Torr O2 at 550 K for **5 min**
 - Valve cracked
 - O2 flowing from O3 conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24, 25 W to ramp (in 20 sec intervals); let climb to 550 K at 25 W (started timer when it hit 550 K); 23.5 W laser power to stabilize at 550 K
 - maintained T within 550-554 K over duration of oxidation
 - Cooling in O2 till T = 310 K (~1 min) , then evacuating HP cell till p_baratron = 0.1 mTorr **10:55 AM**

Temperature profile recorded with QMS: 200203-file1 (T at beginning 285 K)

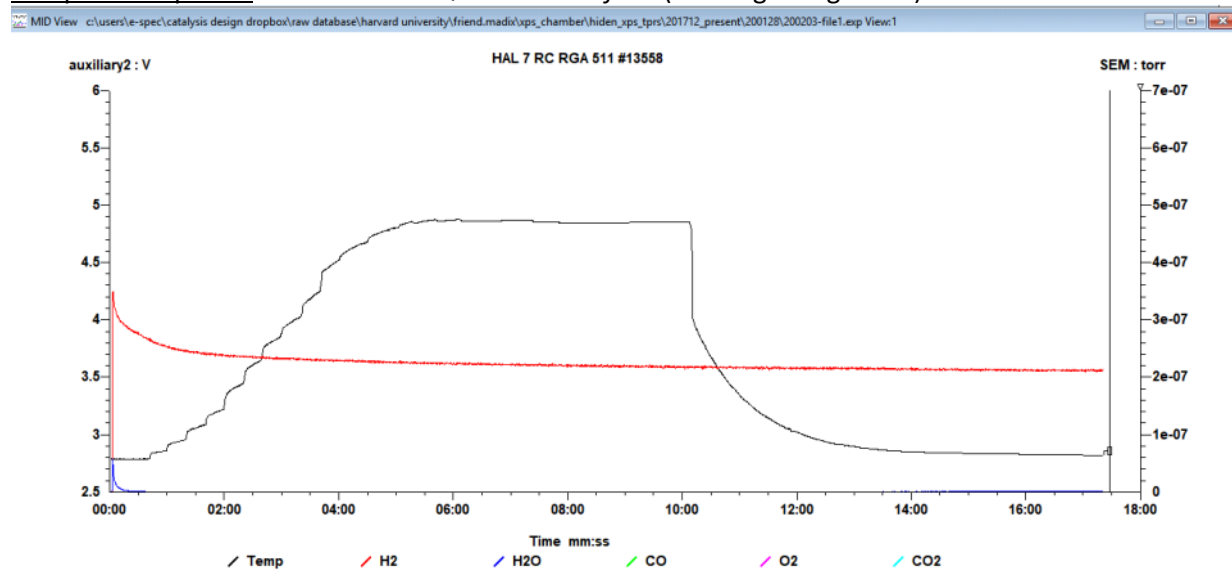


Switch back to air cooling for H2 dose; clear water from cryostat
 Shut valve to HP cell at 15.0"
 Turn on X-ray source at 15.0"




O2 pre-stage pressure was 1925 psi at 11:15 AM

Measure



XPS: after oxidation - 66 degrees


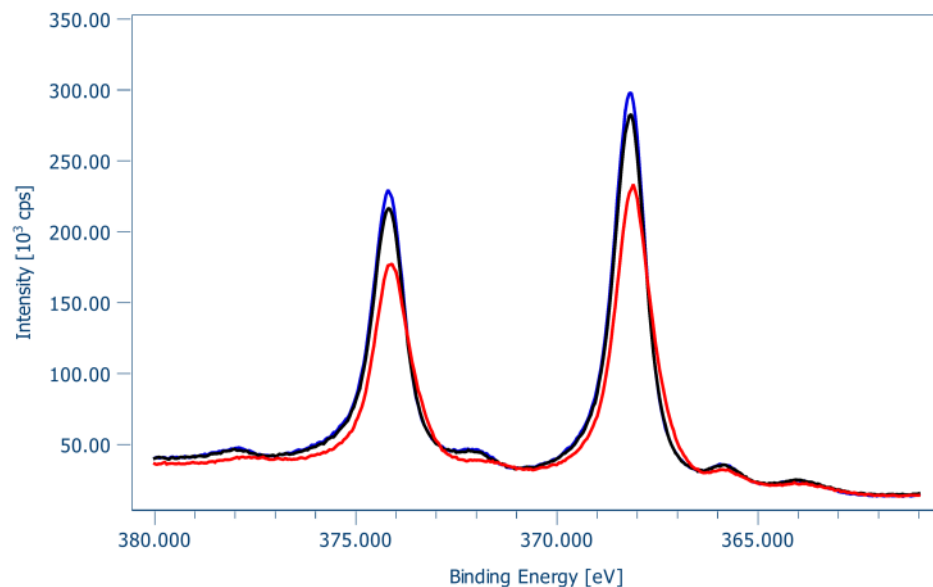


Pd metallic ~335.1-2 eV

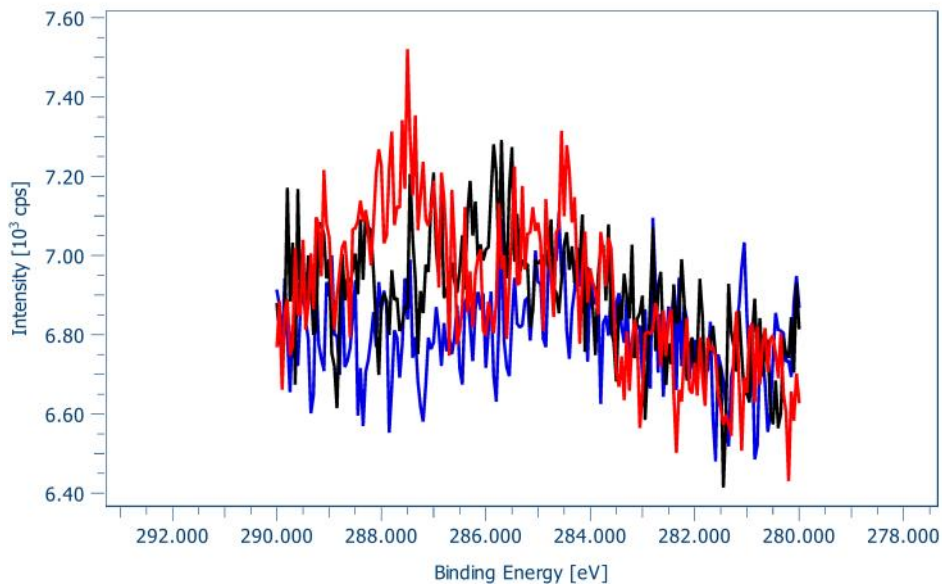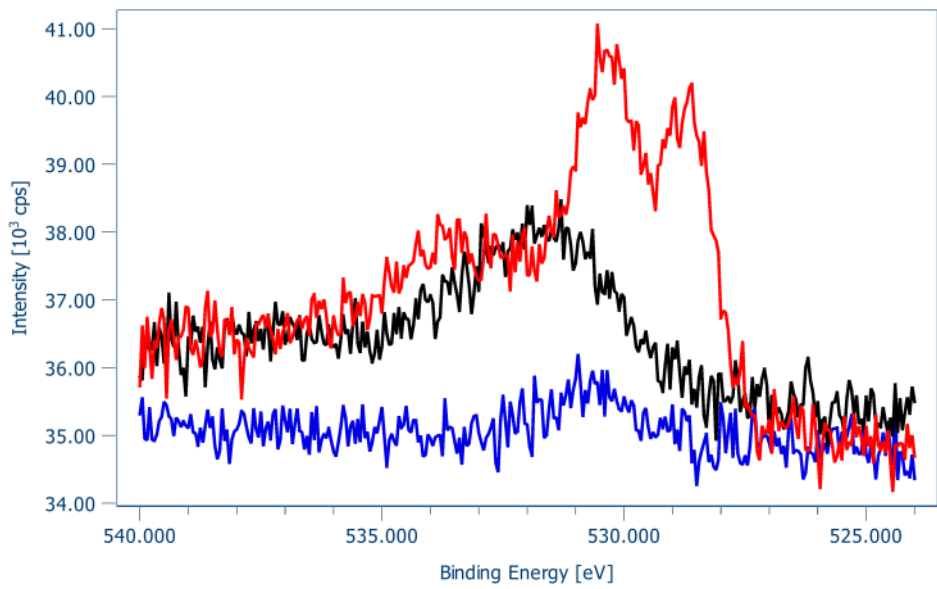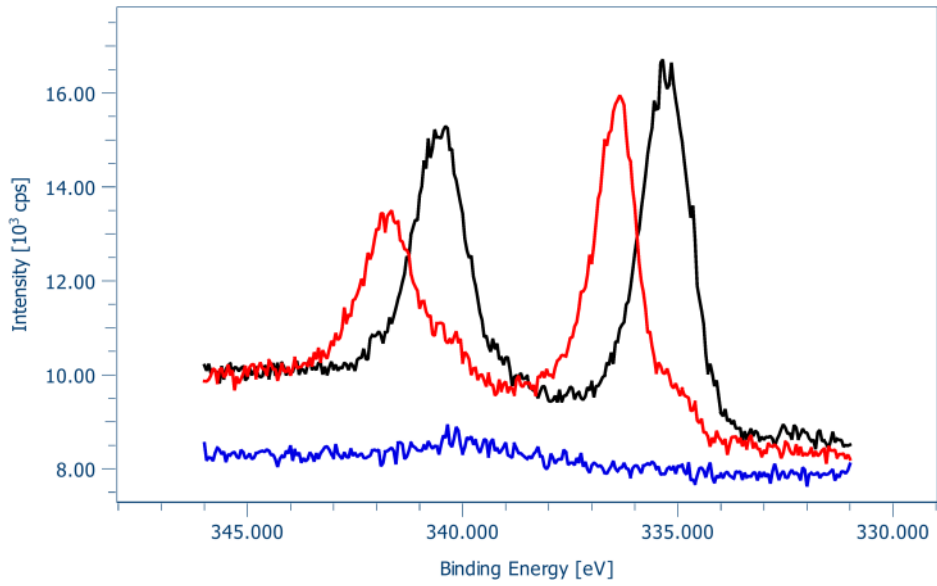

Comparison to 200128 below: Fantastic!, oxidized (red) extremely similar
Reminder 200128 was 0.14 ML Pd!

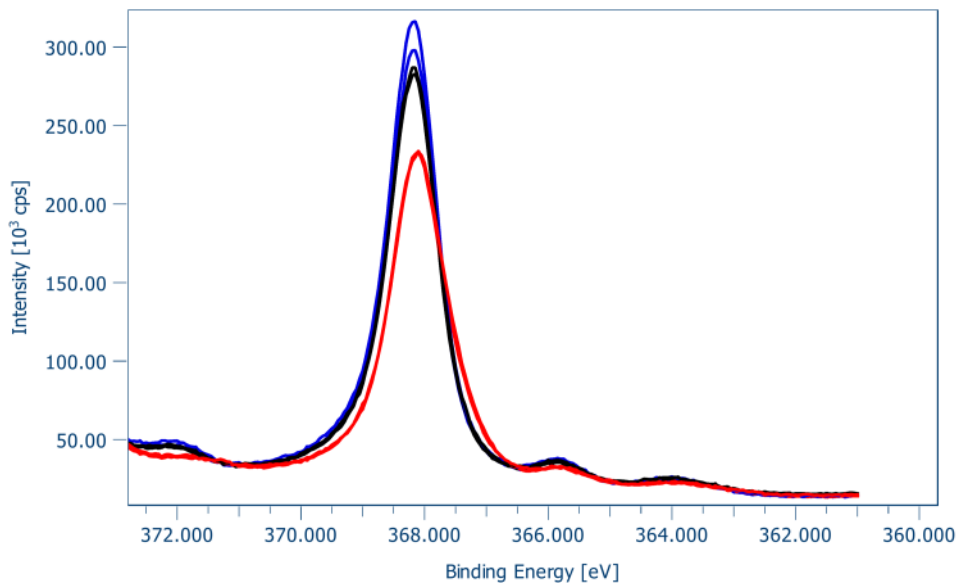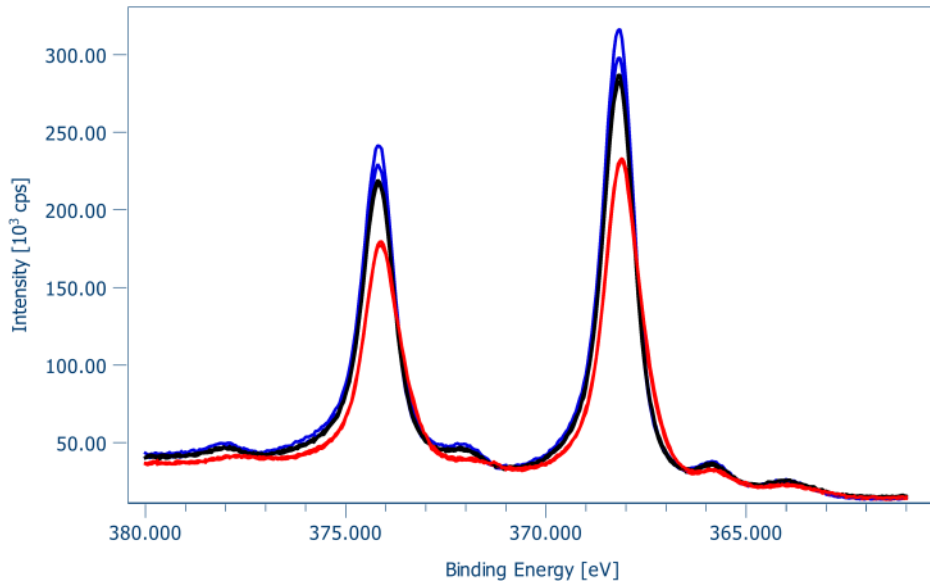

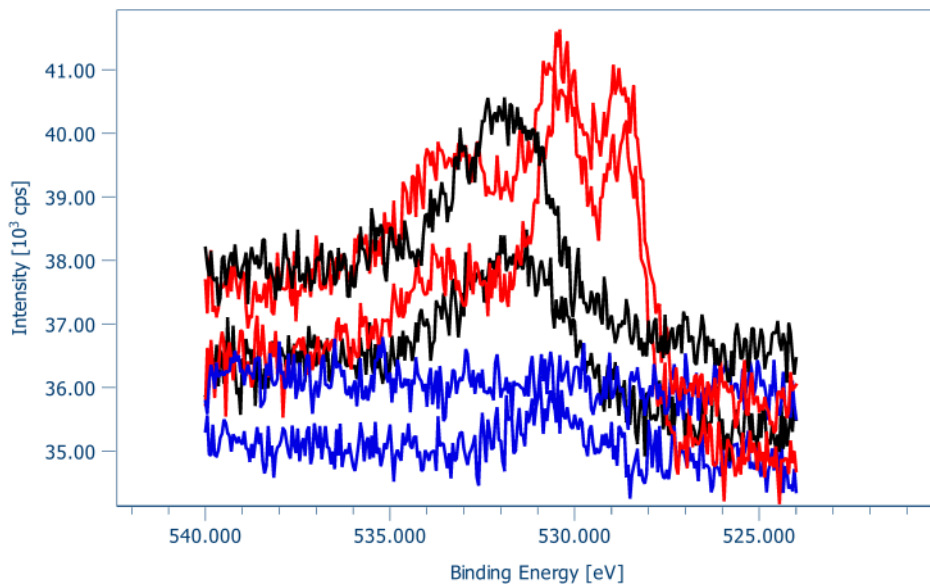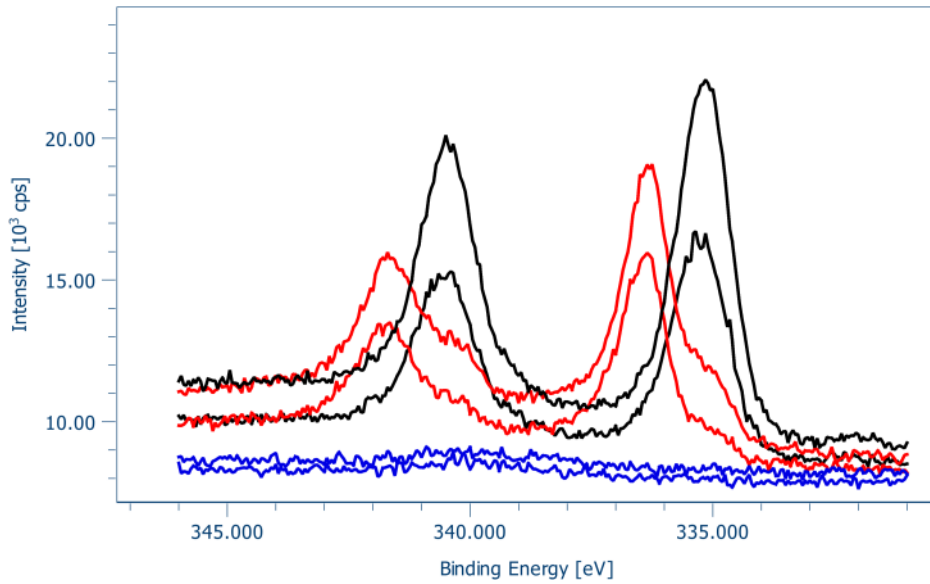

30 mins exposure to background with X-ray gun on

XPS: 30 mins bkg - 66 degrees

Start sputter 12:12 AM

Start O3 12:10 AM

- Sputter for 45 mins at 20 mA and 3e-5 Torr Ar
- Anneal to 800 K for 10 min

#####

2nd round

#####

Leak O3 in HP cell 1:00 PM

Stop leak O3 in HP cell 2:05 PM

XPS: clean Ag#2 - 66 degrees

Switch on evap

Connect cooling

2:20 PM

- **Pd evaporation:** 2.50 A, 1000 V, 10.0 mA, 1.6-1.2 nA for **1 min 30 sec** at 300K
 - Goal: 0.10 ML
 - $0.14/1.5 = 0.0933$
 - $0.1/0.0933 = 1.0718$
 - $1.0718 * 2.0/1.7 = 1.2609$
 - $0.26 * 60 = 15.6$

Switch on X-Ray source before transfer

XPS: Pd deposition#2 - 0 degrees

- Pd area: 10 kcps eV
- Ag area: 717 kcps eV
- Pd/Ag: $10/717 = 0.0139$ --> **0.08(4) ML fine compared to deposition before**

XPS: Pd deposition#2 - 66 degrees


- Pd area: 8.5 kcps eV
- Ag area: 280 kcps eV
- Pd/Ag: $8.5/280 = 0.0304$ --> **0.08(3) ML very accurate angle at low coverages**

Done at 3:01 PM

Turn up O2 at 3:14 APM

- **Oxidized** in 1.92(6) Torr O2 at 550 K for **5 min**
 - Valve cracked
 - O2 flowing from O3 conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24, 25, 25.5 W to ramp (in 20 sec intervals); let climb to 550 K at 25.5 W (started timer when it hit 550 K); 23.5 W laser power to stabilize at 550 K
 - maintained T within 550-554 K over duration of oxidation
 - Cooling in O2 till T = 302 K (~1 min) , then evacuating HP cell till p_baratron = 0.1 mTorr **10:55 AM**

Temperature profile recorded with QMS: 200203-file2 (T at beginning 286 K)


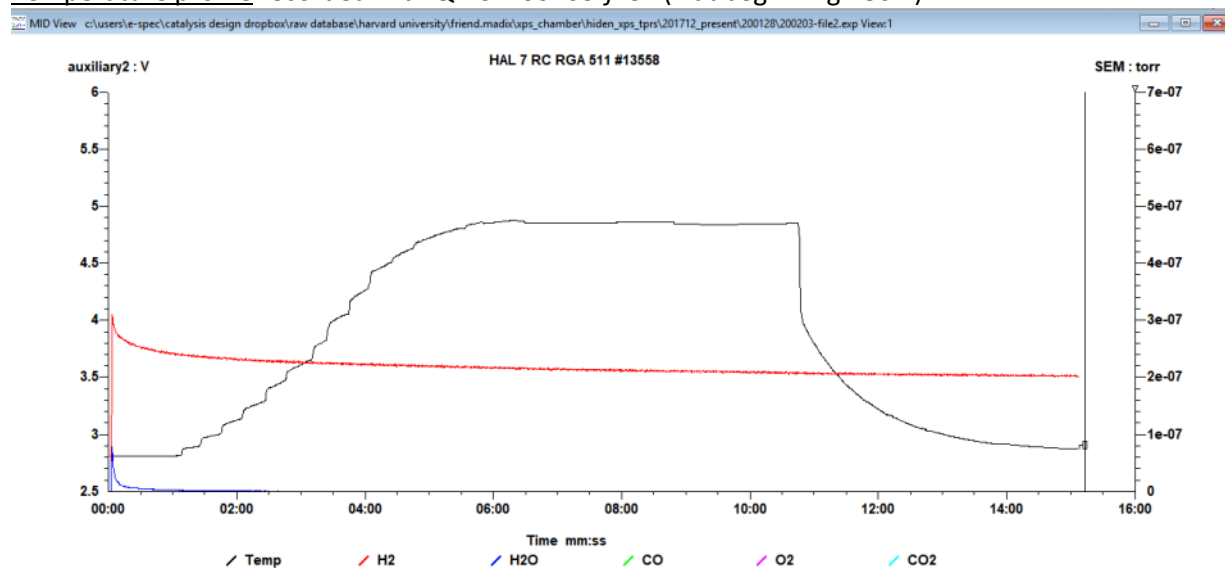




Switch back to air cooling for H2 dose; clear water from cryostat
Shut valve to HP cell at 15.0"
Turn on X-ray source at 15.0"

O2 pre-stage pressure was 1925 psi at 11:15 PM




Measure


XPS: after oxidation#2 - 66 degrees

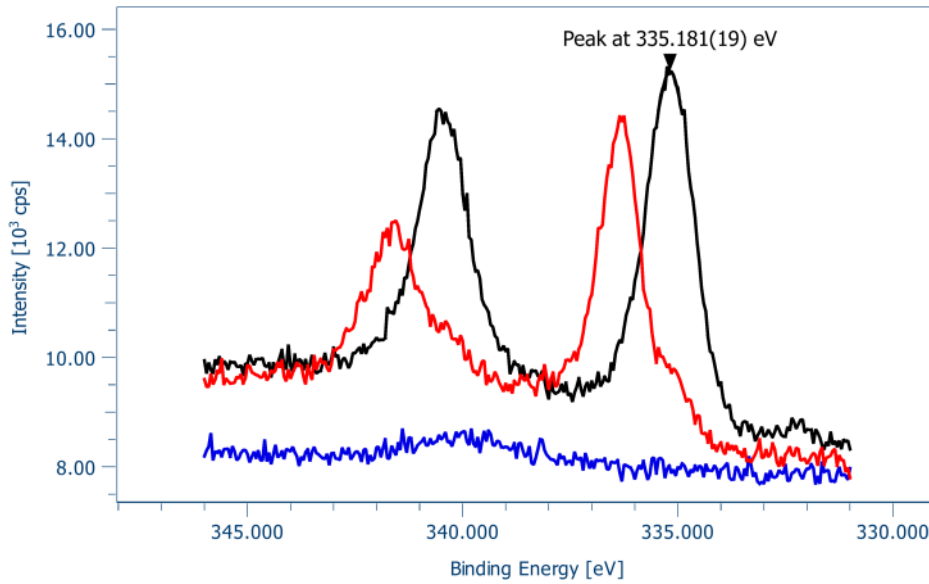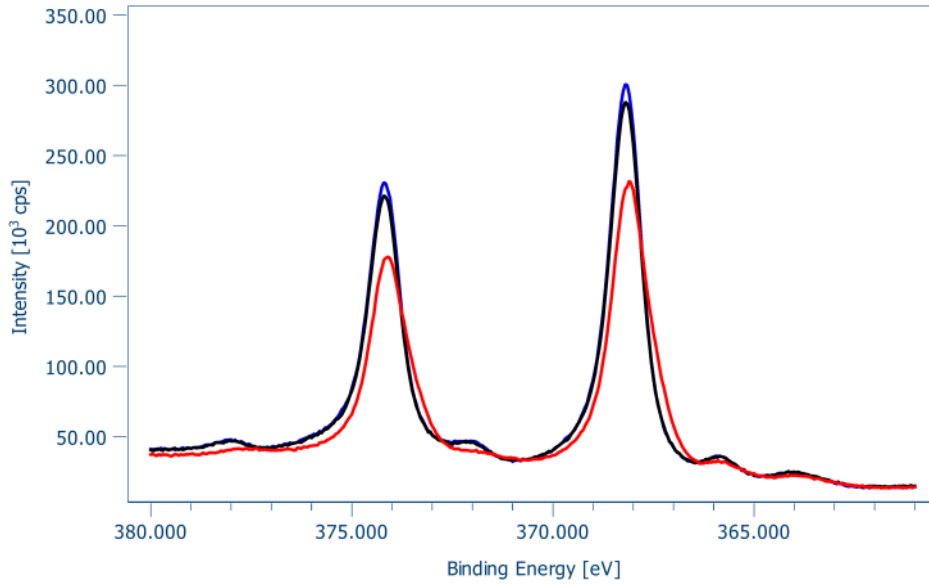

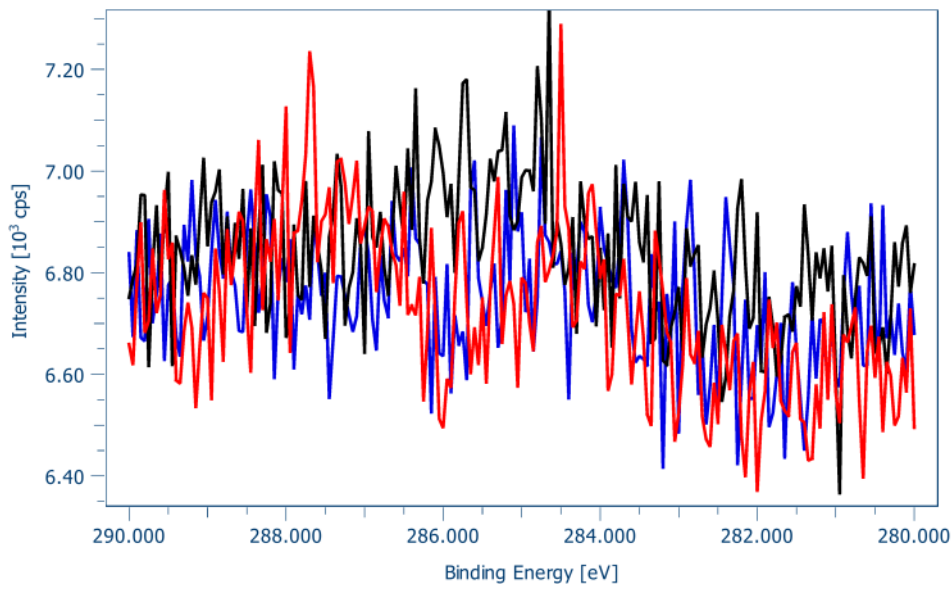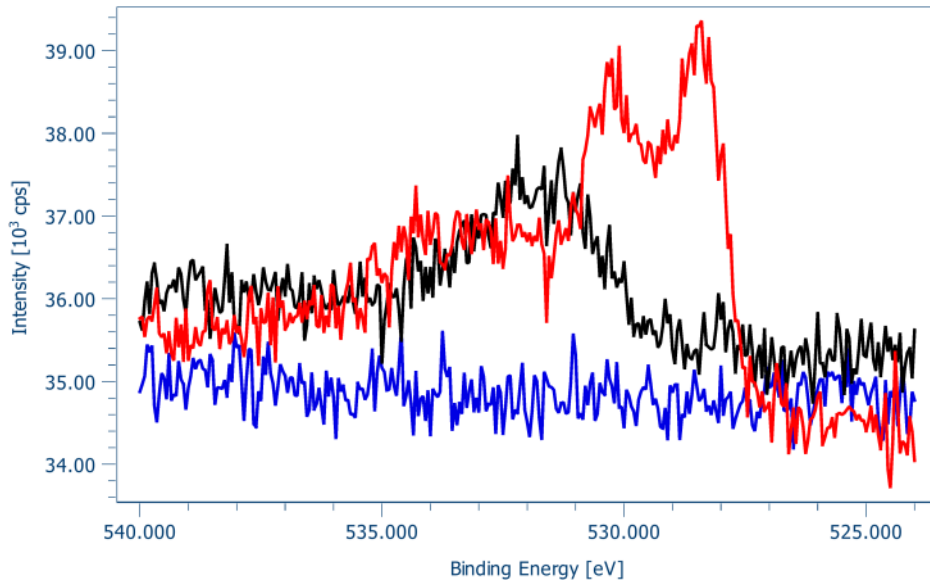


Comparison to first oxidation today: 2nd round in dashed lines
Calling this reproducible

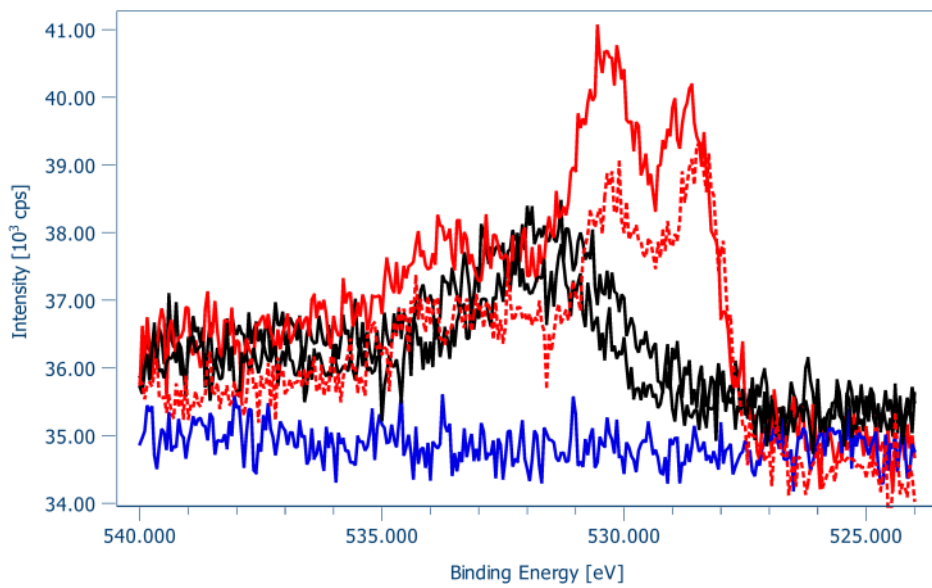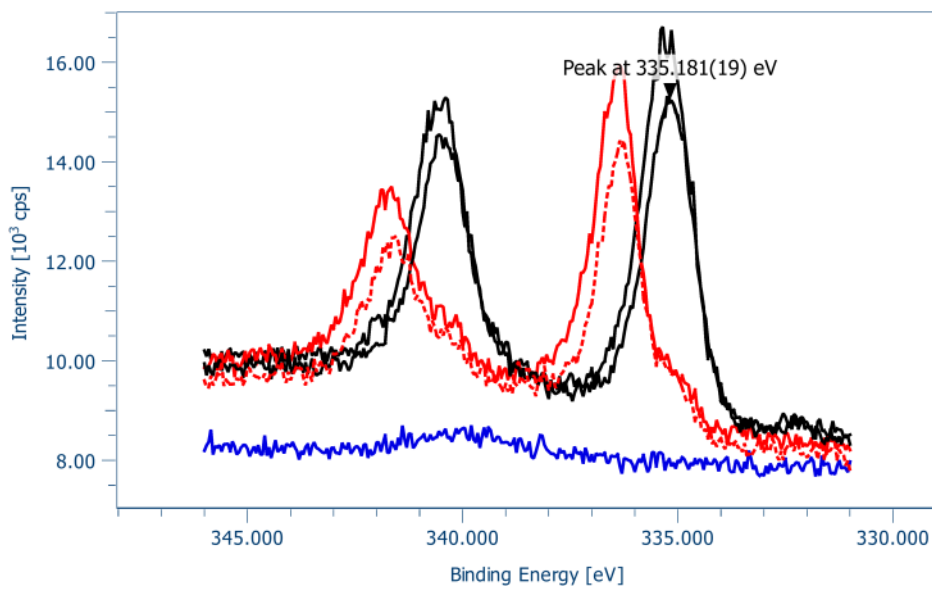
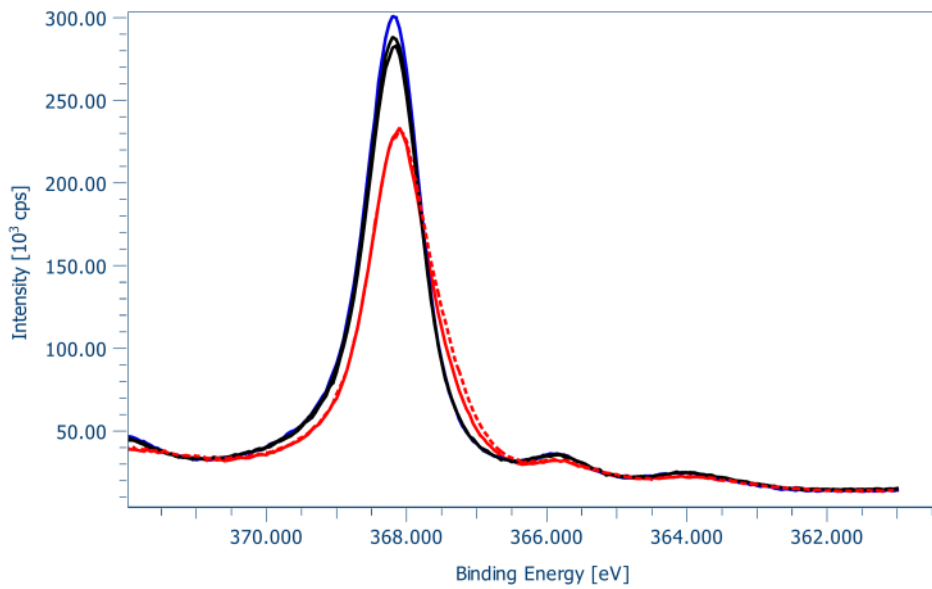


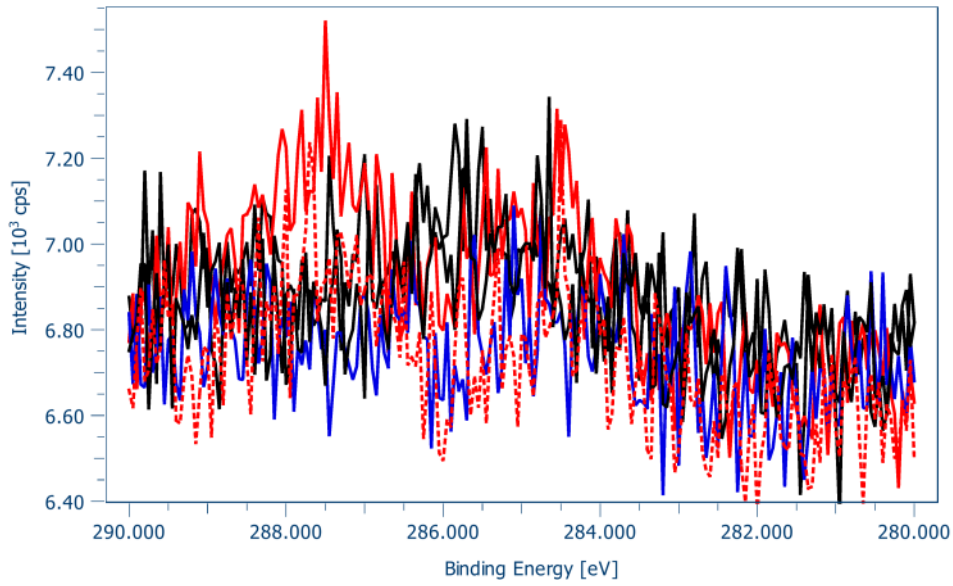

30 mins exposure to background with X-ray gun off

XPS: 30 mins bkg#2 - 66 degrees
purple

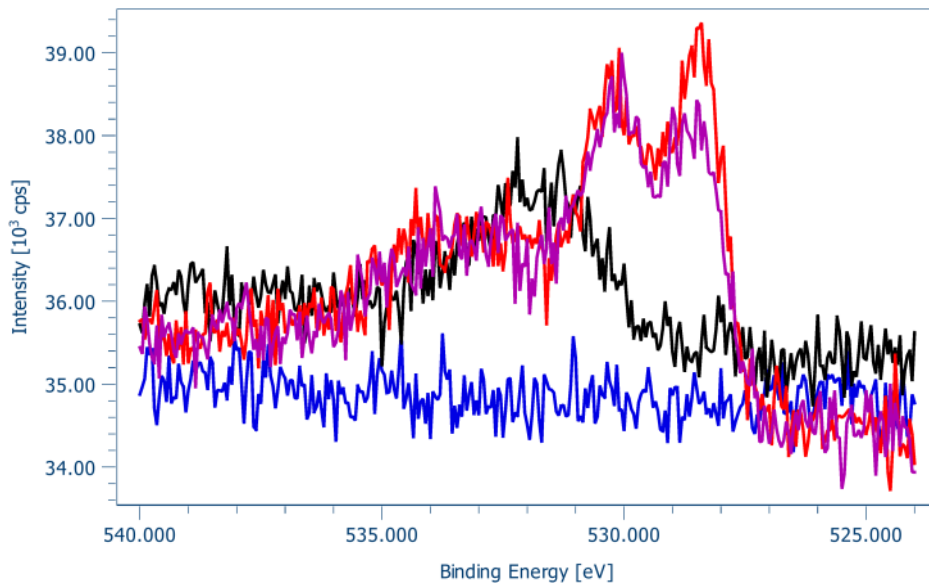

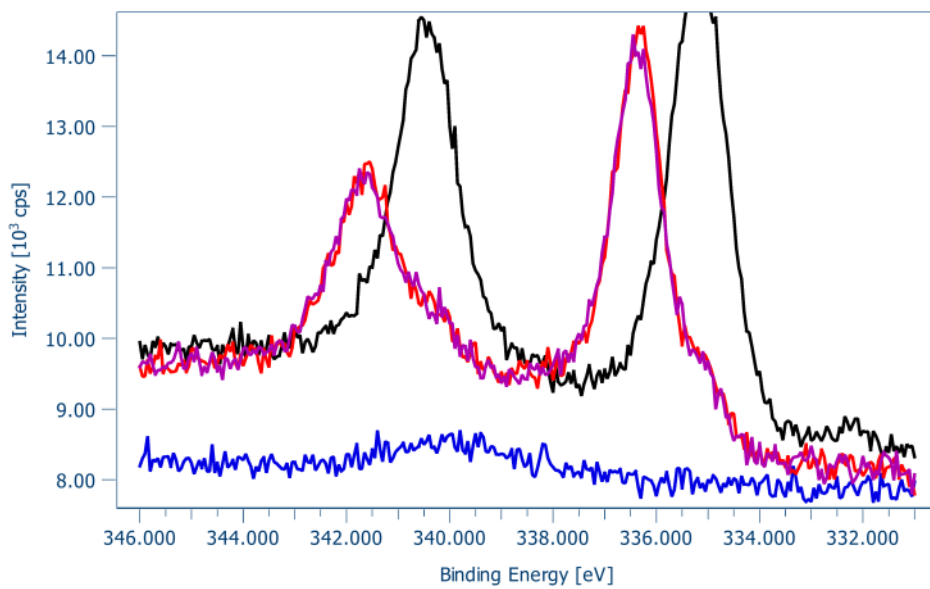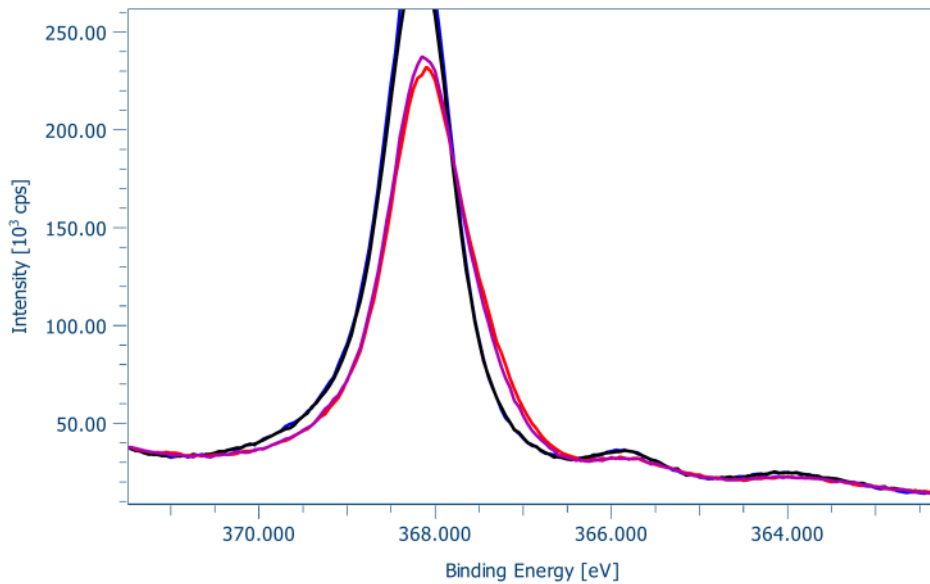

Start sputter 5:00 PM

- Sputter for 30 mins at 20 mA and 3e-5 Torr Ar
- Anneal to 800 K for 10 min

Conclusion:

State of chamber not changing

Established protocol which has to be followed! (#EXPTL SUM / 200203 - O2 Protocol.docx)

0.09 or 0.14 ML doesn't make a difference as expected

200204 - PdAg + 575 & 600 K O2 + H2

Tuesday, February 4, 2020 9:53 AM

Plan

Oxidation temp where PdO 100% oxidized

Logbook

Legend - Plan in purple

Flash anneal to 800 K

Set up water cooling

Open gate valve

Measure clean Ag and turn down O3 and off (8:50 AM)

XPS: clean Ag - 66 degrees

Switch on evap

- **Pd evaporation:** 2.50 A, 1000 V, 9.7 mA, 1.1 nA for **2 min 00 sec** at 300K
Goal: 0.12 ML
 $0.12/0.0533=2.2514$

Connect cooling

9:55 AM

Switch on X-Ray source before transfer and turn ion gauge off

XPS: Pd deposition - 0 degrees

- Pd area: 15 kcps eV
- Ag area: 702 kcps eV
- Pd/Ag: $15/702 = 0.0214$ --> **0.12(8) ML GREAT**

XPS: Pd deposition - 66 degrees


- Pd area: 14 kcps eV
- Ag area: 268 kcps eV
- Pd/Ag: $14/268 = 0.0522$ --> **0.14(0) ML TOP**

Done at 10:30 AM

Turn up O2 at 10:45 AM

- **Oxidized** in 1.93(2) Torr O2 at 575 K for **5 min**
 - Valve cracked
 - O2 flowing from O3 conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 25 W to ramp (in 20 sec intervals); let climb to 575 K at 25 W (started timer when it hit 575 K); 23.5 W laser power to stabilize at 550 K
 - maintained T within 575-579 K over duration of oxidation
 - Cooling in O2 till T = 305 K (~1 min) , then evacuating HP cell till p_baratron = 0.1 mTorr **11:12 AM**

Temperature profile recorded with QMS: 200204-file1 (T at beginning 285 K)

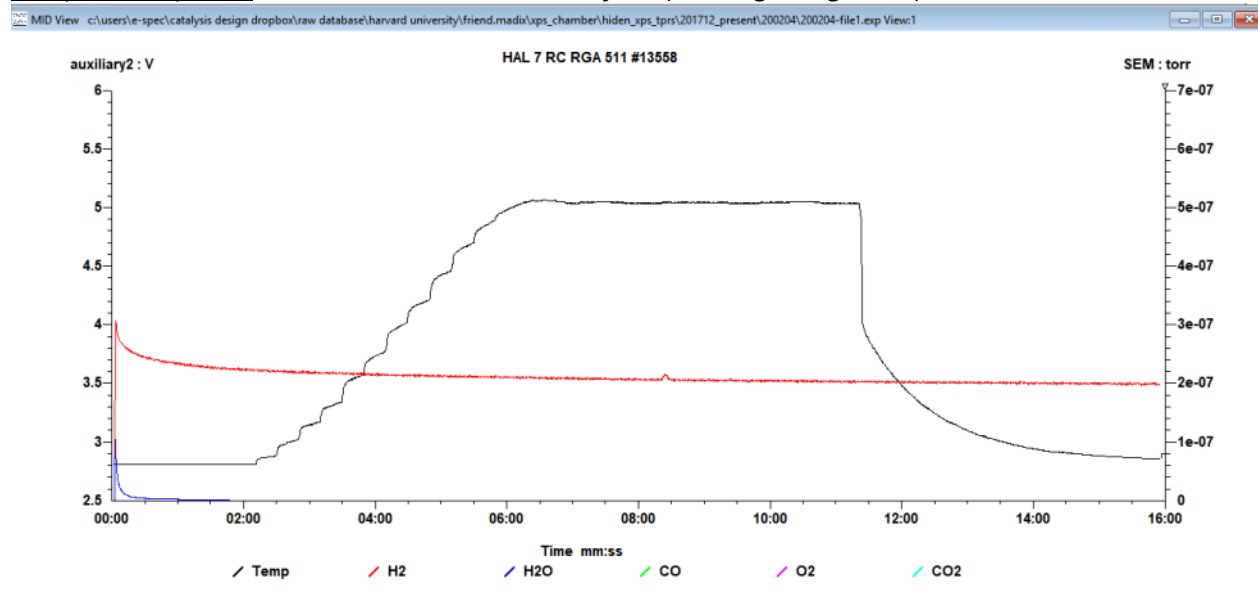

Switch back to air cooling for H2 dose; clear water from cryostat
Shut valve to HP cell at 15.0"
Turn on X-ray source at 15.0"


O2 pre-stage pressure was 1575 psi at 11:12 AM



Measure (11:22 AM)

XPS: after oxidation - 66 degrees

Ag seems a bit less oxidized than 550 K

Pd oxidation pretty similar to 550 K

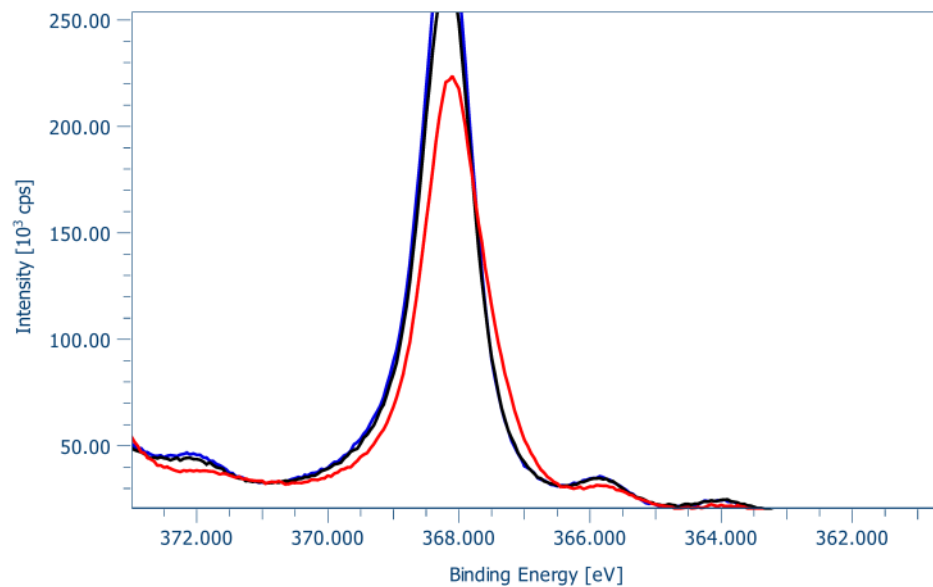

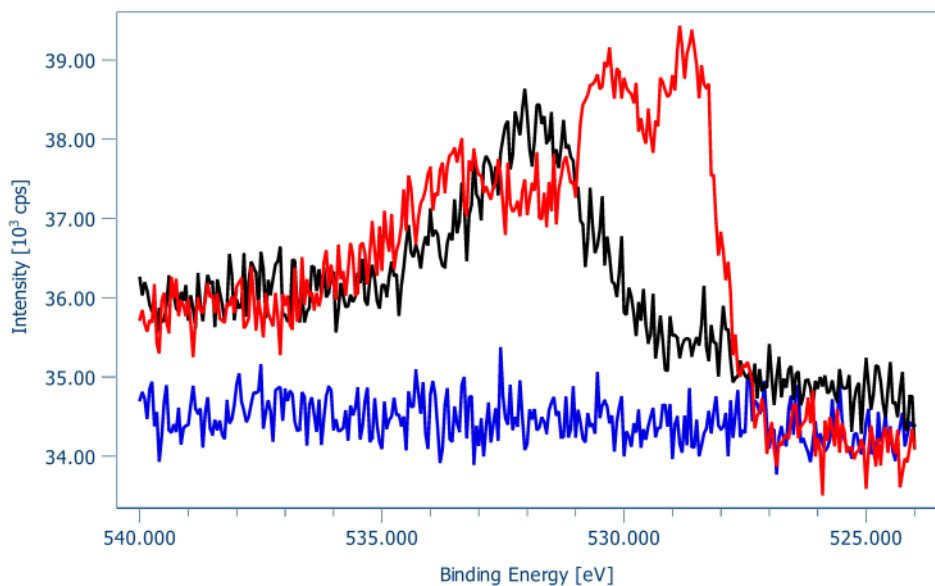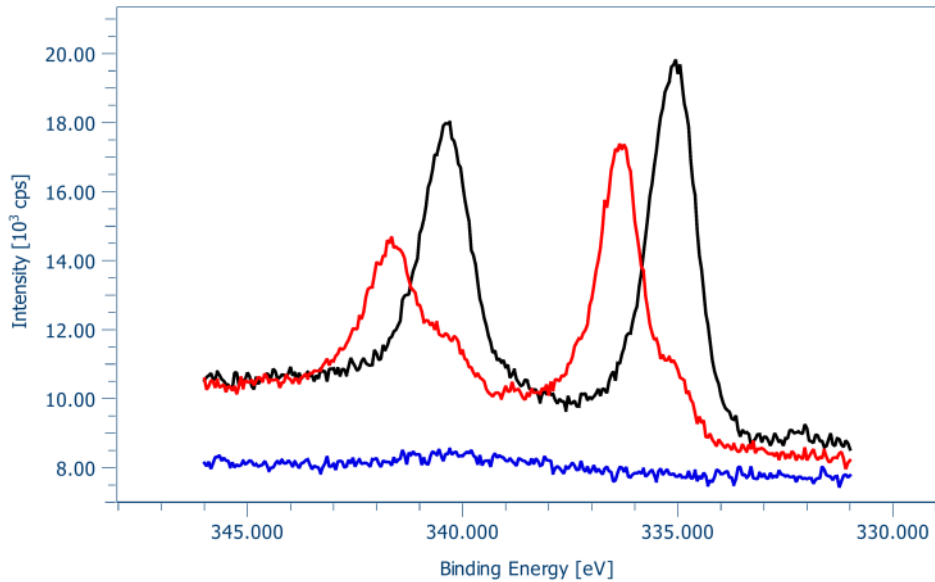

XPS: after oxidation - 0 degrees

Pd buried in deeper layers?

Done at 11:44 AM


2nd oxidation (subsequent)
 #####

Turn up O2 at 12:20 PM

- **Oxidized in 1.93(7) Torr O2 at 600 K for 5 min**
 - Valve cracked
 - O2 flowing from O3 conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 25, 27, 29 W to ramp (in 20 sec intervals); let climb to 600 K at 29 W (started timer when it hit 600 K); 27.8 W laser power to stabilize at 600 K
 - maintained T within 599-604 K over duration of oxidation

- Cooling in O₂ till T = 305 K (~1 min) , then evacuating HP cell till p_baratron = 0.1 mTorr 12:40 PM
- H₂O cooling was very efficient this time

Temperature profile recorded with QMS: 200204-file2 (T at beginning 282 K)

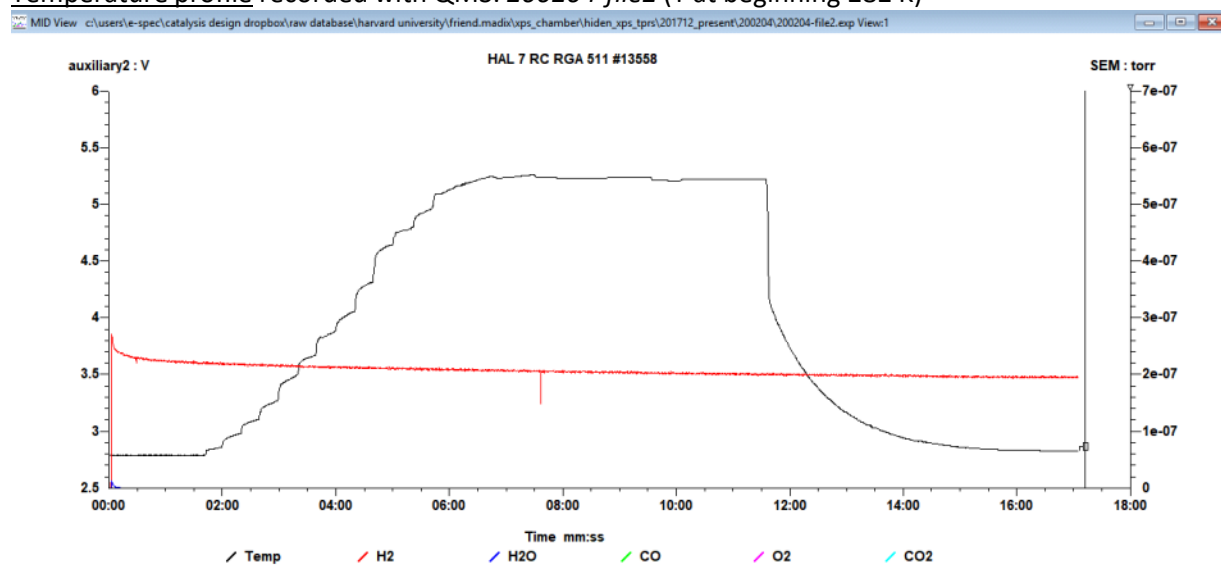


Switch back to air cooling for H₂ dose; clear water from cryostat
 Shut valve to HP cell at 15.0"
 Turn on X-ray source at 15.0"

Measure (12:50 PM)

XPS: after oxidation - 66 degrees

Almost equal compared to previous 550 K oxidation

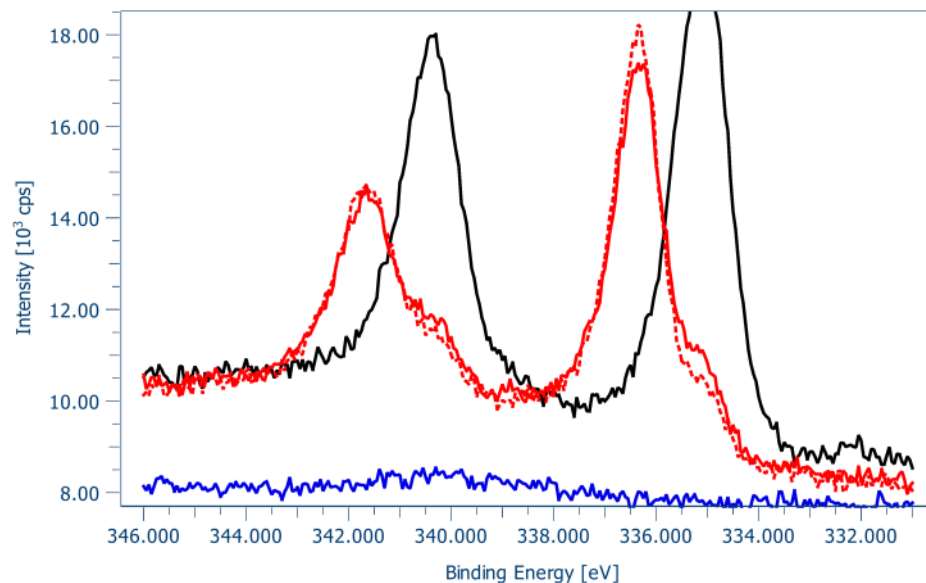

- **H₂ Dose#1:** 2e-4 H₂ for 5 min at 300 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 XPS: H₂ Dose#1 - 66 degrees
- **H₂ Dose#2:** 2e-4 H₂ for 15 min (total: 20 min) at 300 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)

XPS: H2 Dose#2 - 66 degrees

- **H2 Dose#3:** 2e-4 H2 for 15 min (total: 35 min) at 300 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

XPS: H2 Dose#3 - 66 degrees

- **H2 Dose#4:** 2e-4 H2 for 15 min (total: 50 min) at 300 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

XPS: H2 Dose#4 - 66 degrees

Pd deposition



After oxidation#2



H2 Dose#1

H2 Dose#2

H2 Dose#3

H2 Dose#4

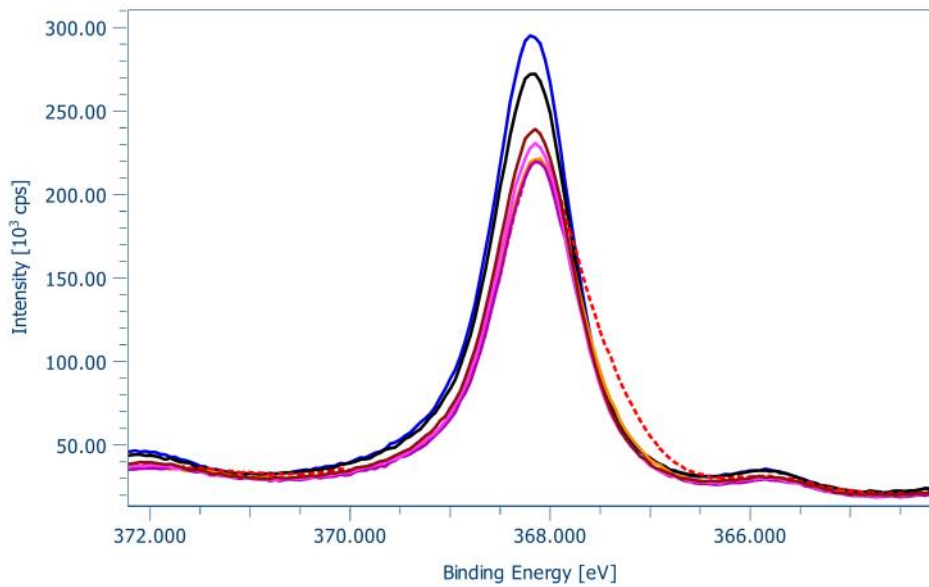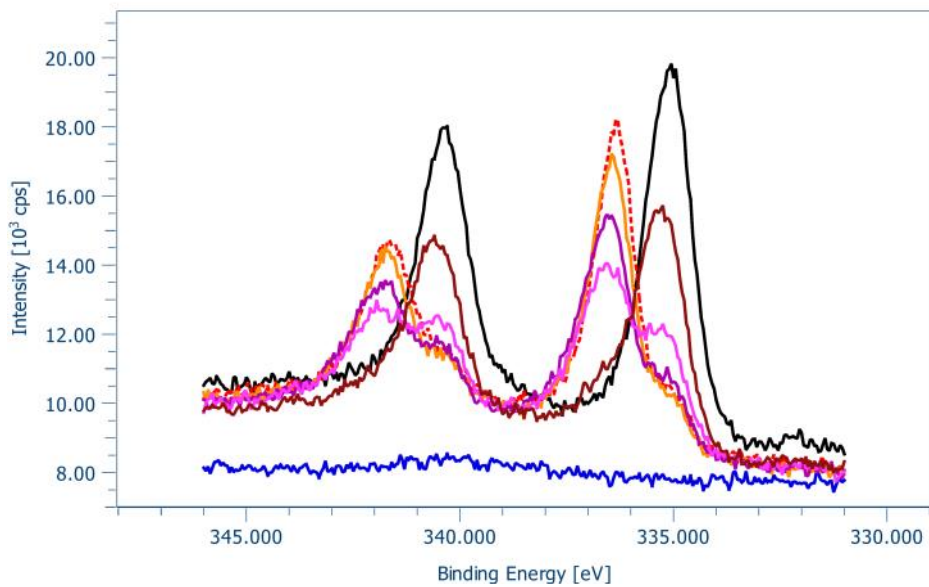

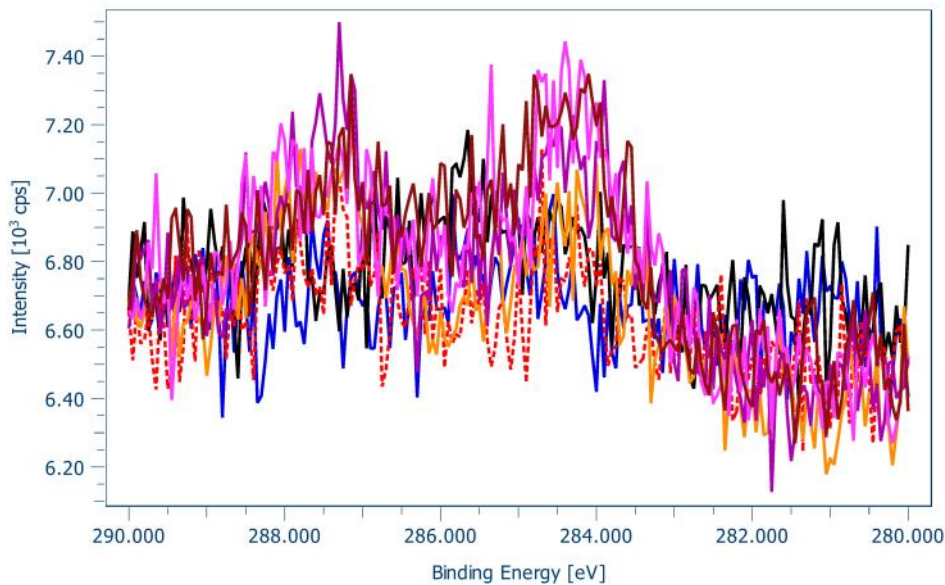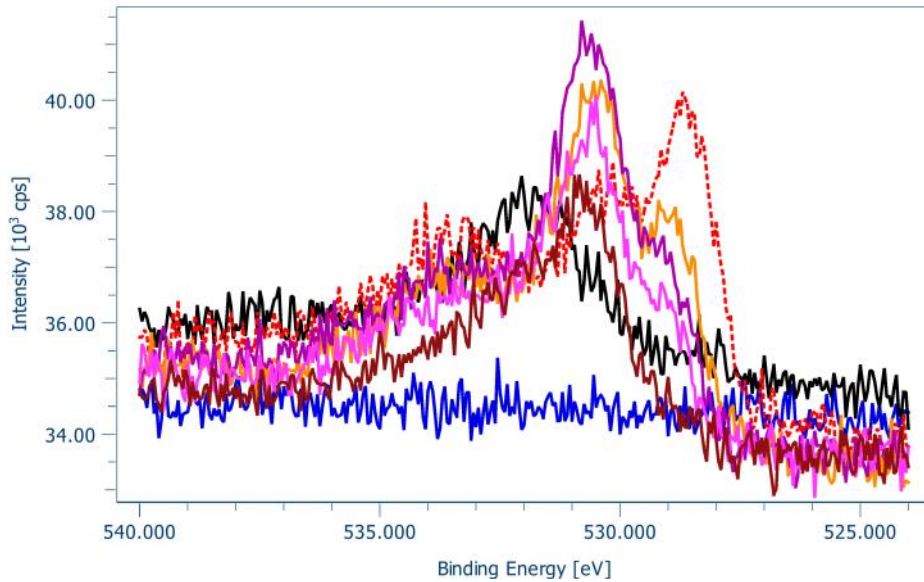

Conclusion

1. Differences of higher (600 K) and lower (550 K) T oxidation are small; choose rather 550 K since this is closer to conditions making p4x4
2. Time steps for Pd at E-4 Torr H2: 5, 20, 35, 40, 45, 50 + min
3. Additional time steps for Ag at E-4 Torr H2: 1, 3 min
4. Do also 1 min at E-7 Torr H2

200206 - PdAg + 550 K O2 + H2 [1st cycle]

Thursday, February 6, 2020 9:01 AM

Flash anneal to 800 K
Set up water cooling
Open gate valve

Turn down O3 and off (9:00 AM)

Measure clean Ag

XPS: clean Ag - 66 degrees

Connect cooling

10:00 AM

- **Pd evaporation:** 2.50 A, 1000 V, 9.6 mA, 1.2 - 0.6 nA for **2 min 00 sec** at 300K
Goal: 0.12 ML

Switch on X-Ray source before transfer and turn ion gauge off

XPS: Pd deposition - 0 degrees

- Pd area: 8.3 kcps eV
- Ag area: 686 kcps eV
- Pd/Ag: $8.3/686 = 0.0121$ --> **0.07(3) ML Not enough**

- **Pd evaporation#2:** 2.50 A, 1000 V, 9.6 mA, 0.7 nA for **1 min 30 sec** at 300K
Goal: 0.12 ML
 $0.07/2 = 0.035$ ML/min
 $0.12 - 0.07 = 0.05$ ML
 $0.05/0.035 = 1.4286$ min
 $0.43 * 60 = 25.8$ sec

XPS: Pd deposition#2 - 0 degrees

- Pd area: 14 kcps eV
- Ag area: 686 kcps eV
- Pd/Ag: $14/686 = 0.0204$ --> **0.12(2) ML GREAT**

Start at 10:34 AM

XPS: Pd deposition#2 - 66 degrees

- Pd area: 12.6 kcps eV
- Ag area: 262 kcps eV
- Pd/Ag: $12.6/262 = 0.0481$ --> **0.12(9) ML NICE**


Done at 10:54 AM

Turn up O2 at 11:04 AM

- **Oxidized** in 1.93(1) Torr O2 at 550 K for **5 min**
 - Valve cracked
 - O2 flowing from O3 conditions before
 - turbo 19 W

- 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24 W to ramp (in 20 sec intervals); let climb to 575 K at 25 W (started timer when it hit 550 K); 20.5 W laser power to stabilize at 550 K
- maintained T within 550-558 K over duration of oxidation
- Cooling in O₂ till T = 305 K (~1 min) , then evacuating HP cell till p_{baratron} = 0.1 mTorr 11:36 AM

Temperature profile recorded with QMS: 200206-file1 (T at beginning 285 K)

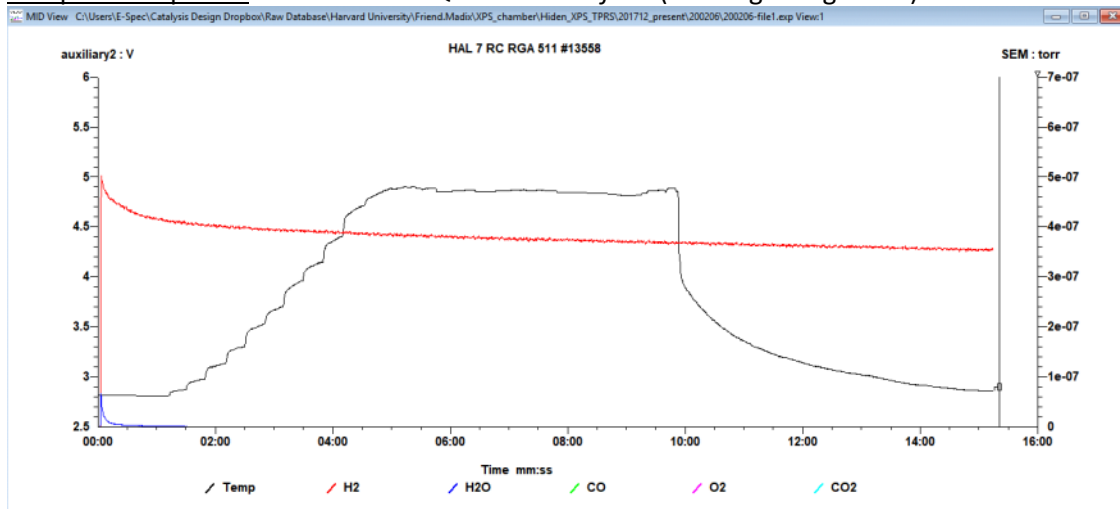


Switch back to air cooling for H₂ dose; clear water from cryostat
 Shut valve to HP cell at 15.0"
 Turn on X-ray source at 15.0"

O₂ pre-stage pressure was 825 psi at 11:36 AM

Measure (11:47 AM)

XPS: after oxidation - 66 degrees

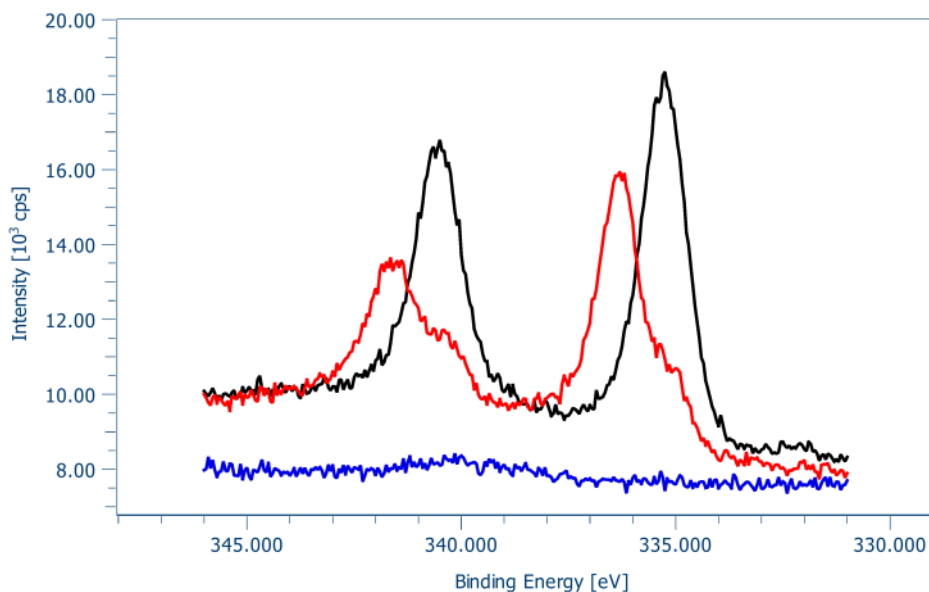

Closed O₂ valves, stop O₂ flow behind HP O₂ leak valve

- **H₂ Dose#1:** 1e-7 H₂ for 1 min at 295 K (background, right leak valve)
 (5e-8 in prep; had analysis ion gauge off)
XPS: H₂ Dose#1 - 66 degrees

p_prep = 3.7e-9
Done 12:09 PM

- **H2 Dose#2:** 1e-5 H2 for 1 min at 296 K (background, right leak valve)
(5.0e-6 in prep; had analysis ion gauge off)
XPS: H2 Dose#2 - 66 degrees
p_prep = 6e-9 - 4.7e-9
Done 12:22 PM
- **H2 Dose#3:** 1e-4 H2 for 1 min at 297 K (background, right leak valve)
(5.0e-5 in prep; had analysis ion gauge off)
Wait 2.5 min to recover pressure
XPS: H2 Dose#3 - 66 degrees
p_prep = 1.9e-8
Done 12:39 PM
- **H2 Dose#4:** 2e-4 H2 for 1 min (total: 1 min) at 298 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 2.5 min to recover pressure
XPS: H2 Dose#4 - 66 degrees
p_prep = 2.6e-8
Done 12:53 PM
- **H2 Dose#5:** 2e-4 H2 for 1 min (total: 2 min) at 300 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 2.5 min to recover pressure
XPS: H2 Dose#5 - 66 degrees
p_prep = 2.0e-8
Done 1:11 PM
- **H2 Dose#6:** 2e-4 H2 for 3 min (total: 5 min) at 300 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#6 - 66 degrees
p_prep = 2.5e-8
Done 1:29 PM
- **H2 Dose#7:** 2e-4 H2 for 3 min (total: 8 min) at 300 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#7 - 66 degrees
p_prep = 2.1e-8
Done at 1:49 PM
Regime of autocatalytic reduction of PdO
- **H2 Dose#8:** 2e-4 H2 for 1 min (total: 9 min) at 300 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#8 - 66 degrees
p_prep = 2.1e-8
Done at 2:08 PM
Regime of autocatalytic reduction of PdO

- **H2 Dose#9:** 2e-4 H2 for 1 min (total: 10 min) at 300 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#9 - 66 degrees
 p_prep = 2.2e-8
 Done at 2:25 PM
 Regime of autocatalytic reduction of PdO
- **H2 Dose#10:** 2e-4 H2 for 5 min (total: 15 min) at 300 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#10 - 66 degrees
 p_prep = 2.3e-8
 Done at 2:48 PM
- **H2 Dose#11:** 2e-4 H2 for 5 min (total: 20 min) at 300 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#11 - 66 degrees
 p_prep = 2.3e-8
 Done at 3:11 PM
 Seems like reduction is complete, but doing one more spectra to be sure
- **H2 Dose#12:** 2e-4 H2 for 10 min (total: 30 min) at 300 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#12 - 66 degrees
 p_prep = 2.3e-8
 Done at 3:40 PM
 Overlapping with last spectra (#11) --> DONE

Pd deposition#2

After oxidation

H2 Dose#1

H2 Dose#2

H2 Dose#3

H2 Dose#4

H2 Dose#5

H2 Dose#6

H2 Dose#7




H2 Dose#8


H2 Dose#9

H2 Dose#10

H2 Dose#11

H2 Dose#12

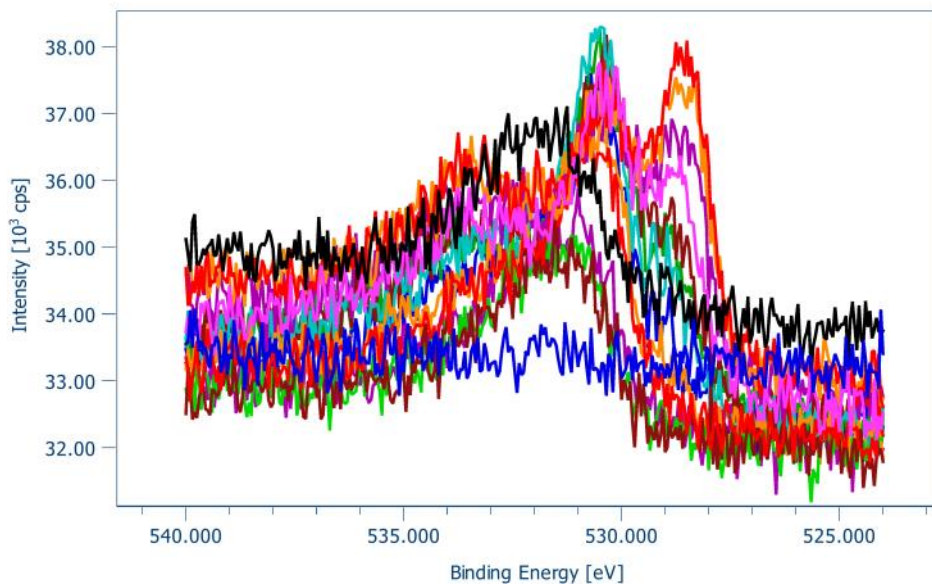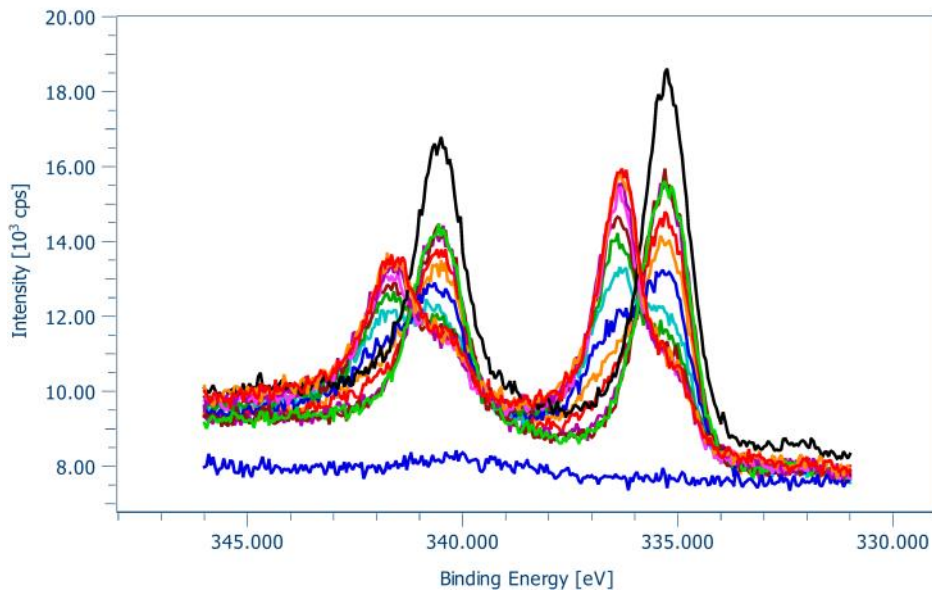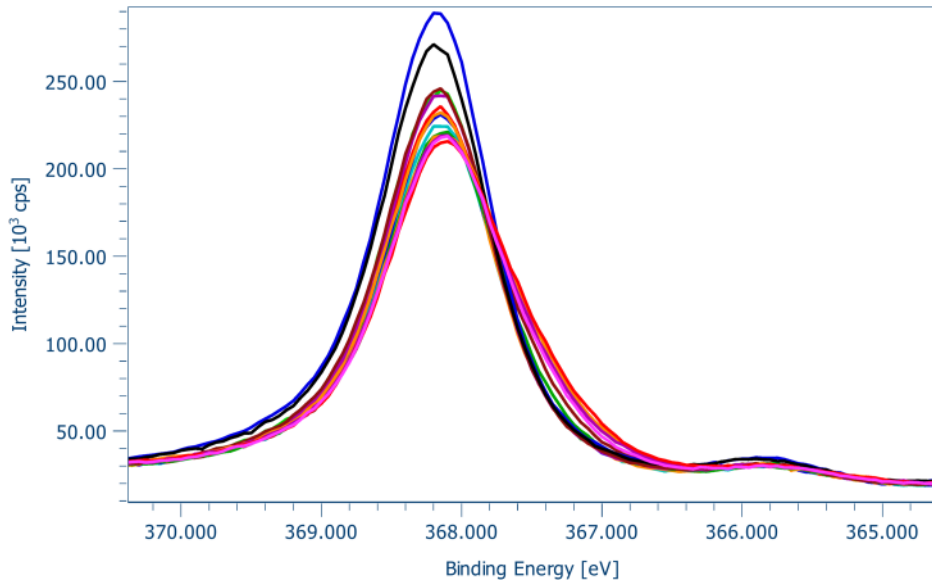

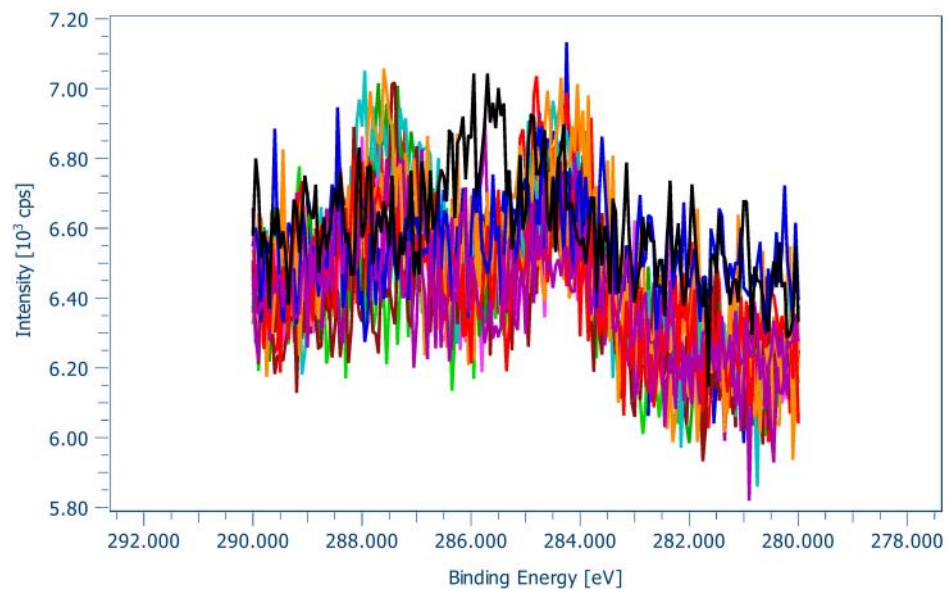

Start O2 flow at 6 PM
Turn up ozone at 6:20 PM

200207 - PdAg + 550 K O₂ + H₂ [2nd cycle]

Friday, February 7, 2020 8:52 AM

Turn down O₃ and off (8:45 AM)

Measure

XPS: Ag after 1st reduction - 66 degrees

Need to leave 1 h 30 mins in UHV (at least till 10:15 AM) between turning down O₃ and oxidizing


Connect cooling

9:10 AM

Turn up O₂ at 10:30 AM

- **Oxidized in 1.92(9) Torr O₂ at 550 K for 5 min**
 - Valve cracked
 - O₂ flowing from O₃ conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24 W to ramp (in 20 sec intervals); let climb to 550 K at 24 W (started timer when it hit 550 K); 21.5 W laser power to stabilize at 550 K
 - maintained T within 550-556 K over duration of oxidation
 - Cooling in O₂ till T = 305 K (~1 min) , then evacuating HP cell till p_{baratron} = 0.1 mTorr 11:00 AM

Temperature profile recorded with QMS: 200207-file1 (T at beginning 284 K)

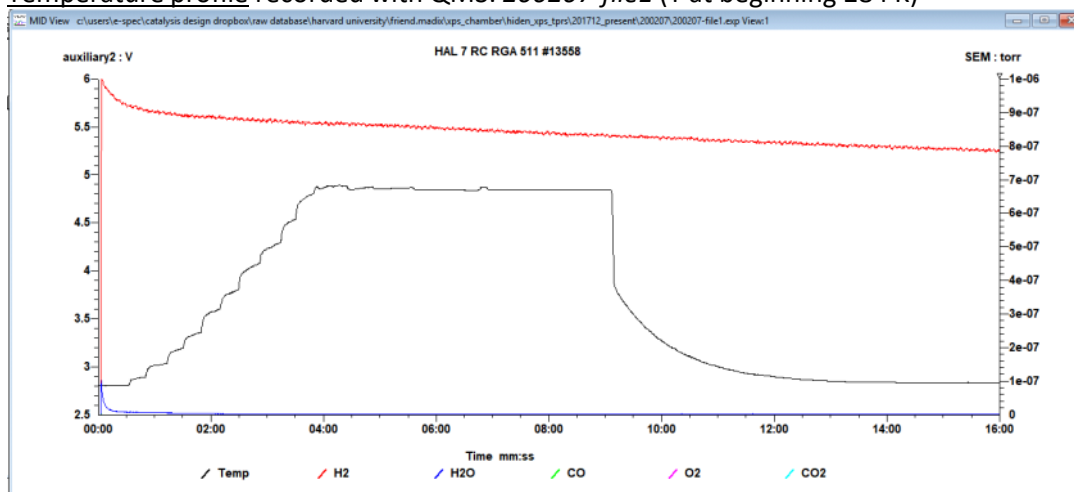

Switch back to air cooling for H₂ dose; clear water from cryostat

Shut valve to HP cell at 15.0"

Turn on X-ray source at 15.0"

O₂ pre-stage pressure was 550 psi at 11:15 AM



Measure (11:11 AM)

XPS: after oxidation - 66 degrees

H₂ Dose#12

Ag after 1st cycle

After oxidation

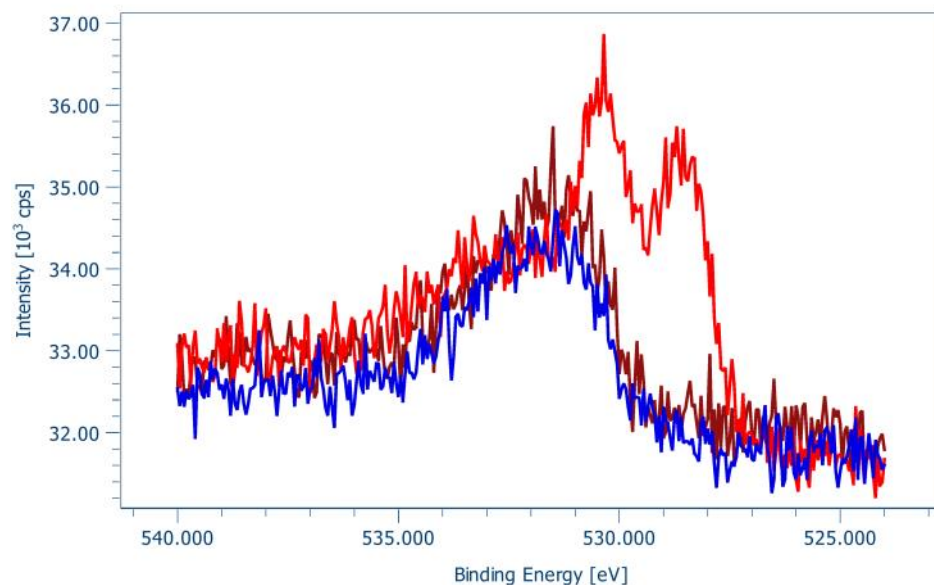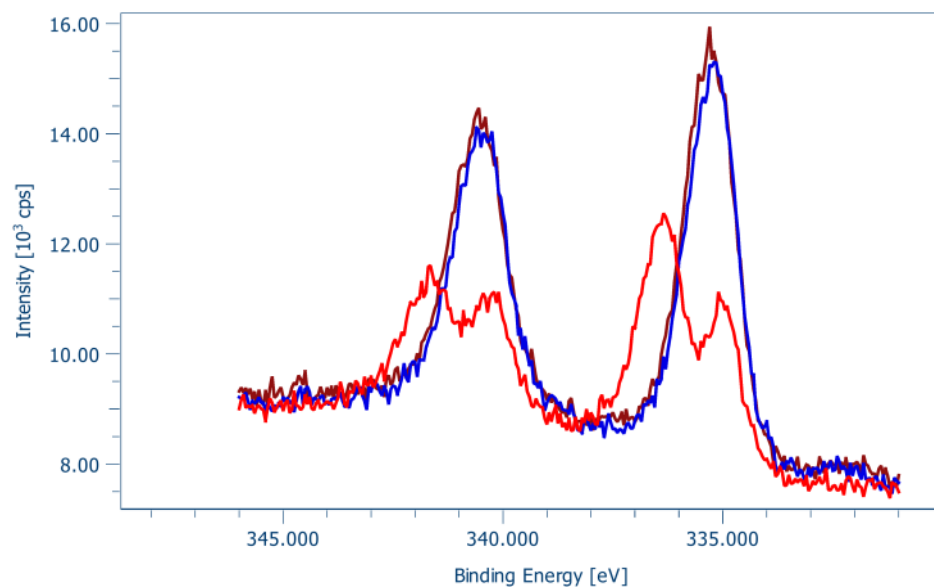

SAME DOSES AS 1ST CYCLE

- **H2 Dose#1:** $1e-7$ H2 for 1 min at 294 K (background, right leak valve)
($5.0e-8$ in prep; $p_{\text{bkg}} = 6e-9$; $p_{\text{delta}} = 4.4e-8$; had analysis ion gauge off)
XPS: H2 Dose#1 - 66 degrees
 $p_{\text{prep}} = 6.3e-9$
Done 11:36 AM
- **H2 Dose#2:** $1e-5$ H2 for 1 min at 295 K (background, right leak valve)
($5.0e-6$ in prep; had analysis ion gauge off)
XPS: H2 Dose#2 - 66 degrees
 $p_{\text{prep}} = 8e-9$
Done 11:48 AM
- **H2 Dose#3:** $1e-4$ H2 for 1 min at 296 K (background, right leak valve)

(5.0e-5 in prep; had analysis ion gauge off)

Wait 2.5 min to recover pressure

XPS: H2 Dose#3 - 66 degrees

p_prep = 2.0e-8

Done 12:06 PM

- **H2 Dose#4:** 2e-4 H2 for 1 min (total: 1 min) at 296 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 2.5 min to recover pressure
XPS: H2 Dose#4 - 66 degrees
p_prep = 2.3e-8
Done 12:22 PM
- **H2 Dose#5:** 2e-4 H2 for 1 min (total: 2 min) at 296 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 2.5 min to recover pressure
XPS: H2 Dose#5 - 66 degrees
p_prep = 2.3e-8
Done 12:38 PM
- **H2 Dose#6:** 2e-4 H2 for 3 min (total: 5 min) at 297 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#6 - 66 degrees
p_prep = 2.2e-8
Done 12:59 PM
- **H2 Dose#7:** 2e-4 H2 for 3 min (total: 8 min) at 297 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#7 - 66 degrees
p_prep = 2.3e-8
Done at 1:19 PM
Regime of autocatalytic reduction of PdO
- **H2 Dose#8:** 2e-4 H2 for 1 min (total: 9 min) at 297 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#8 - 66 degrees
p_prep = 2.1e-8
Done at 1:38 PM
Regime of autocatalytic reduction of PdO
- **H2 Dose#9:** 2e-4 H2 for 1 min (total: 10 min) at 297 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#9 - 66 degrees
p_prep = 2.2e-8
Done at 1:55 PM
Regime of autocatalytic reduction of PdO
- **H2 Dose#10:** 2e-4 H2 for 5 min (total: 15 min) at 297 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#10 - 66 degrees

p_prep = 2.1e-8

Done at 2:19 PM

- **H2 Dose#11:** 2e-4 H2 for 5 min (total: 20 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#11 - 66 degrees

p_prep = 2.5e-8

Done at 2:40 PM

- **H2 Dose#12:** 2e-4 H2 for 10 min (total: 30 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#12 - 66 degrees

p_prep = 2.5e-8

Done at 3:09 PM

- **H2 Dose#13:** 2e-4 H2 for 10 min (total: 40 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Overshot to 2.2e-4 in prep

Wait 3 min to recover pressure

XPS: H2 Dose#13 - 66 degrees

p_prep = 2.6e-8

Done at 3:40 PM

Spectra overlapping with #12, but small bump in O1s --> Do one last dose to make sure

- **H2 Dose#14:** 2e-4 H2 for 10 min (total: 50 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#14 - 66 degrees

p_prep = 2.e-8

Done at 4:13 PM

Overlapping --> DONE!

Don't forget VB spectrum!!

Legend:

H2 Dose#12 (day before)

After oxidation

H2 Dose#1

H2 Dose#2

H2 Dose#3

H2 Dose#4

H2 Dose#5

H2 Dose#6

H2 Dose#7

H2 Dose#8

H2 Dose#9



H2 Dose#10



H2 Dose#11

H2 Dose#12

H2 Dose#13

H2 Dose#14 (not shown)

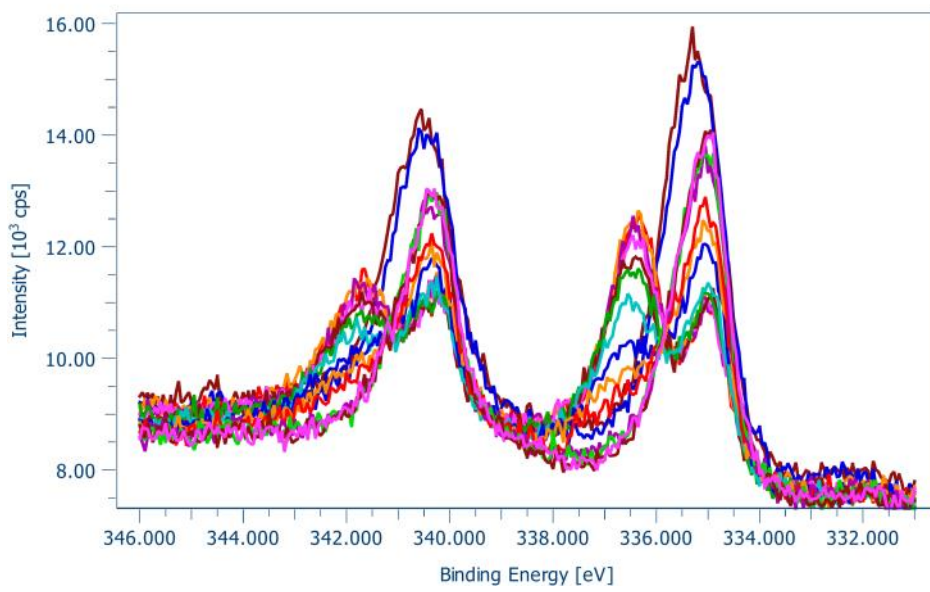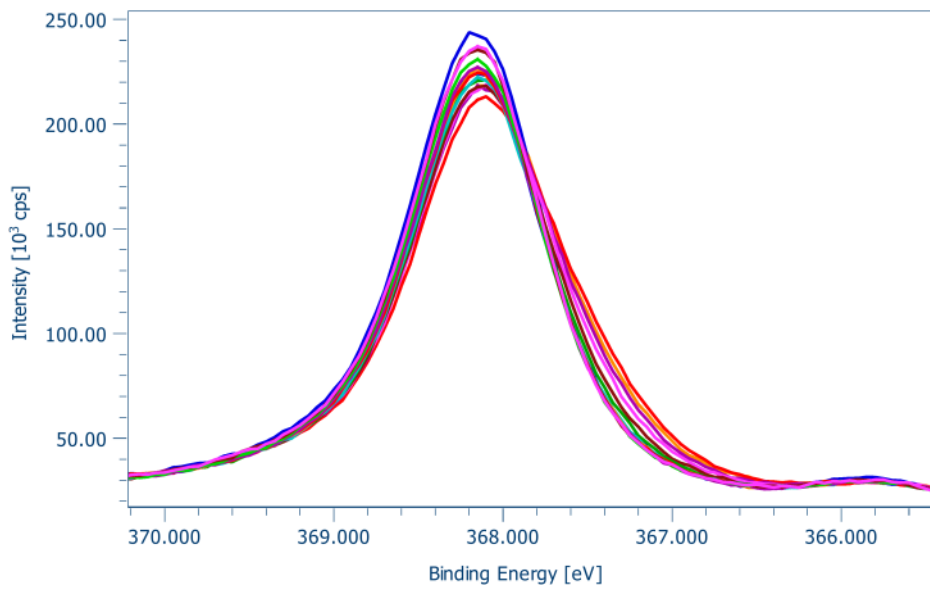

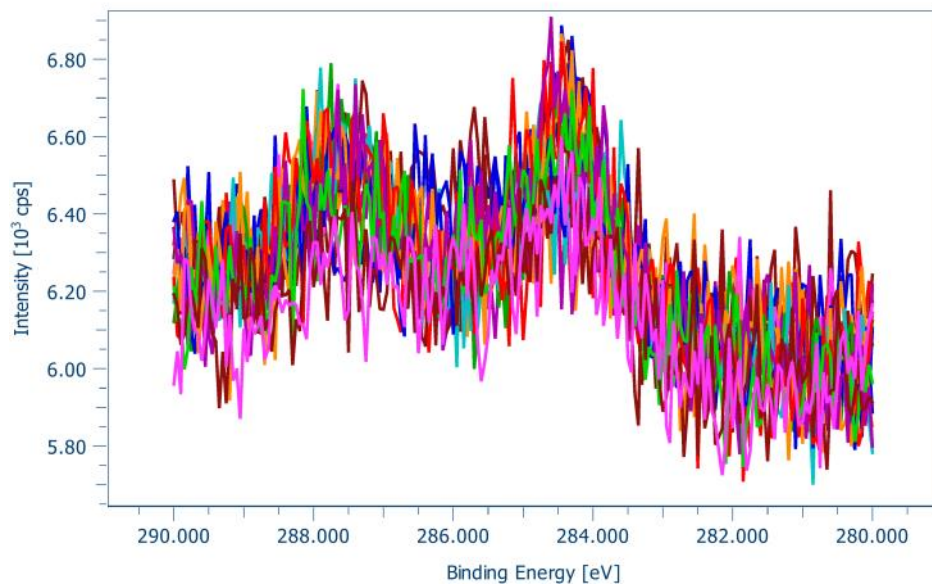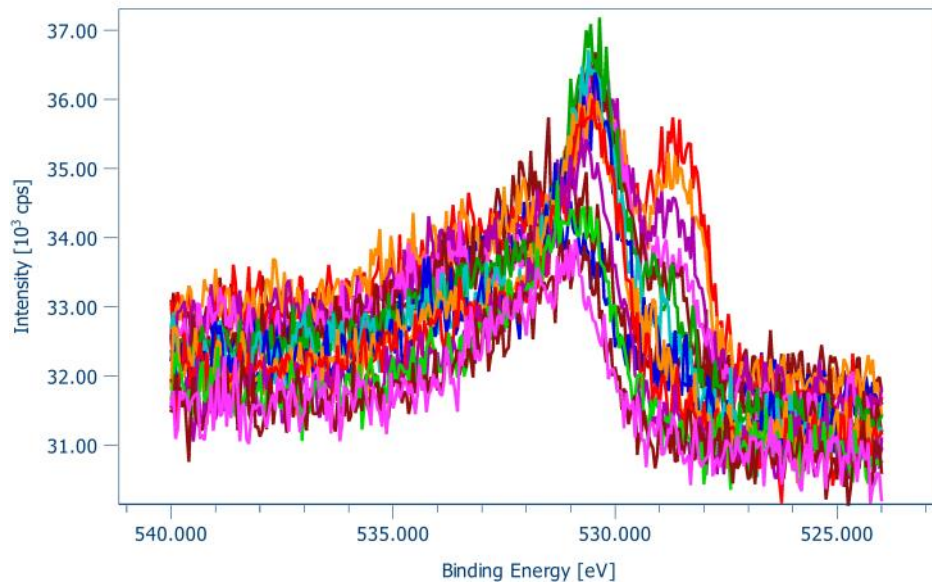

Connected new O2 tank
Pressure in first stage is 2700 psi at 4:00 PM
Turn up O3 at 4:10 PM
Ozone into HP cell 5:45 PM ~30mTorr

200208 - PdAg + 550 K O2 + H2 [3rd cycle]

Saturday, February 8, 2020 9:14 AM

Turn down O3 and off (9:15 AM)

Measure

XPS: Ag after 2nd reduction - 66 degrees
Pressure was $\sim 7e-9$

Need to leave 1 h 30 mins in UHV (at least till 10:45 AM) between turning down O3 and oxidizing


Connect cooling

10:30 AM

Turn up O2 at 11:00 AM

- **Oxidized** in 1.92(6) Torr O2 at 550 K for 5 min
 - Valve cracked
 - O2 flowing from O3 conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24 W to ramp (in 20 sec intervals); let climb to 550 K at 24 W (started timer when it hit 550 K); 21.5 W laser power to stabilize at 550 K
 - maintained T within 550-555 K over duration of oxidation
 - Cooling in O2 till T = 305 K (~ 1 min) , then evacuating HP cell till p_baratron = 0.1 mTorr 11:31 AM

Temperature profile recorded with QMS: 200208-file1 (T at beginning 287 K)

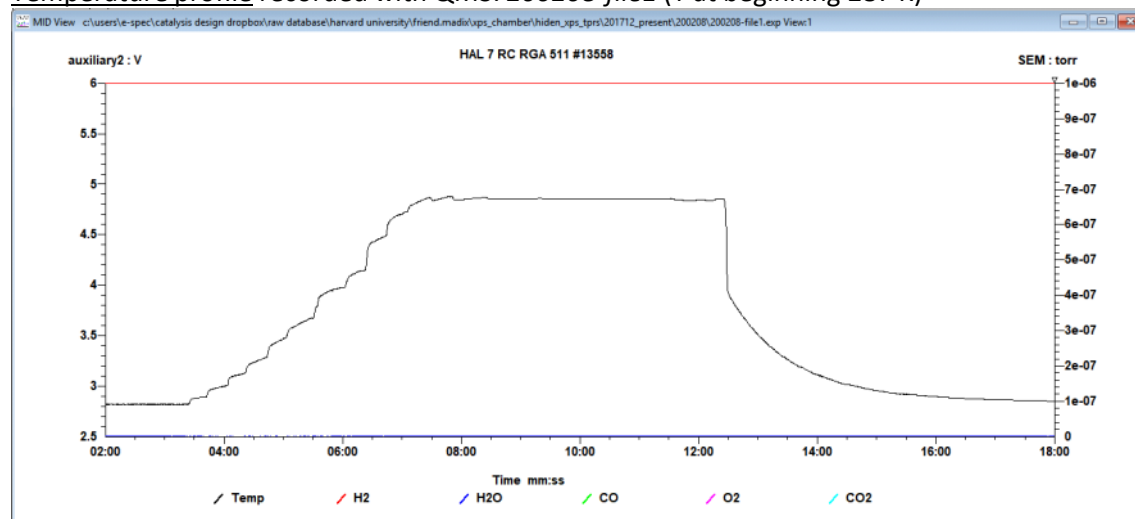

Switch back to air cooling for H2 dose; clear water from cryostat

Shut valve to HP cell at 15.0"

Turn on X-ray source at 15.0"



O2 pre-stage pressure was 2350 psi at 11:44 AM (lost 350 psi)

Closed O2 valves

Measure (11:39 AM)

XPS: after oxidation - 66 degrees

H2 Dose#14
Ag after 2nd reduction
After oxidation

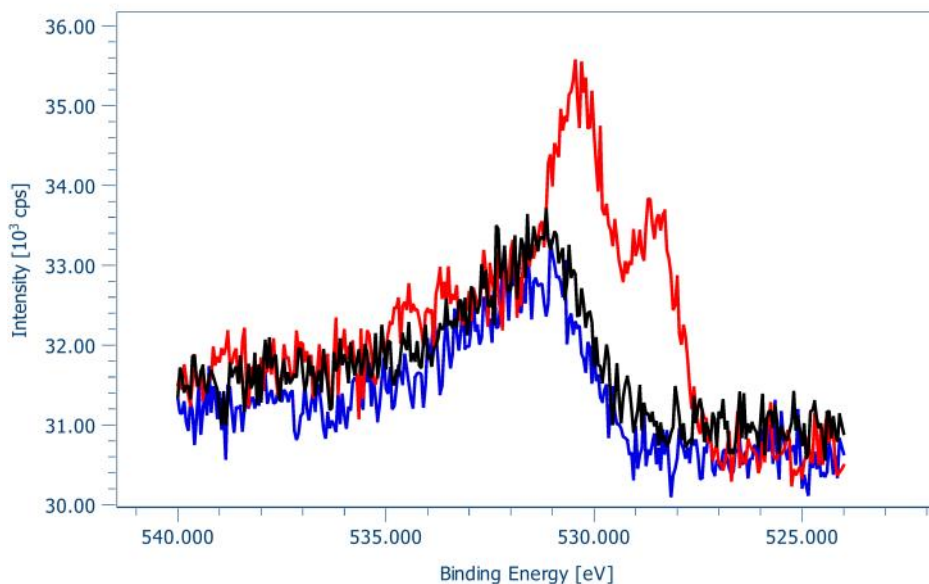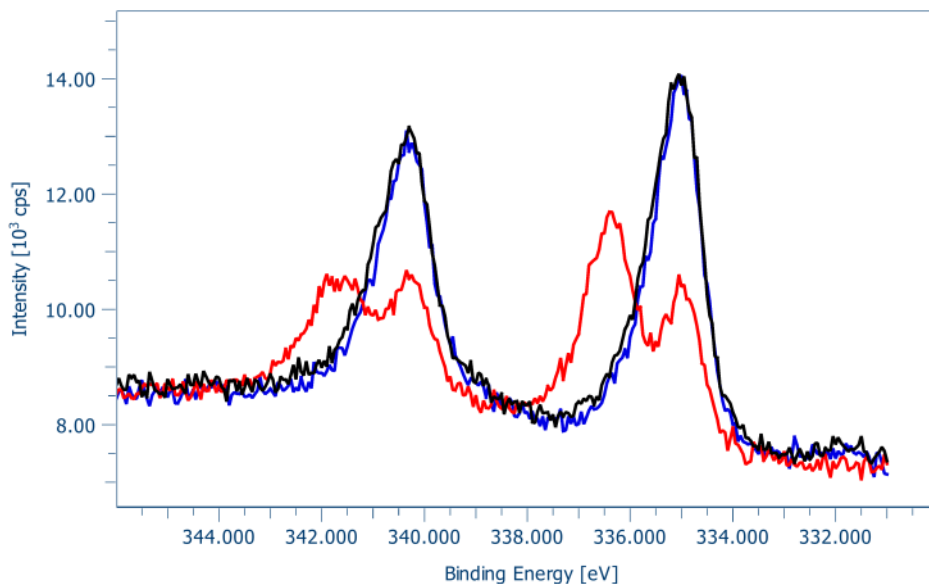

SAME DOSES AS 1ST CYCLE

- **H2 Dose#1:** 1e-7 H2 for 1 min at 295 K (background, right leak valve)
(5.2e-8 in prep; p_{bkg} = 7e-9 p_{delta} = 4.5e-8 ; had analysis ion gauge off)
XPS: H2 Dose#1 - 66 degrees
p_{prep} = 8.3e-9
Done 12:05 PM
- **H2 Dose#2:** 1e-5 H2 for 1 min at 296 K (background, right leak valve)
(5.0e-6 in prep; had analysis ion gauge off)
XPS: H2 Dose#2 - 66 degrees
p_{prep} = 1.0e-8
Done 12:19 PM

- H2 Dose#3:** $1e-4$ H2 for 1 min at 296 K (background, right leak valve)
 (5.0e-5 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure
 XPS: H2 Dose#3 - 66 degrees
 p_prep = 2.1e-8
 Done 12:36 PM
- H2 Dose#4:** $2e-4$ H2 for 1 min (total: 1 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure
 XPS: H2 Dose#4 - 66 degrees
 p_prep = 2.1e-8
 Done 12:54 PM
- H2 Dose#5:** $2e-4$ H2 for 1.5 min (total: 2.5 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure
 Accidentally set timer wrong, did 1.5 min instead of 1 min. Balance out in next dose!
 XPS: H2 Dose#5 - 66 degrees
 p_prep = 2.3e-8
 Done 1:11 PM
- H2 Dose#6:** $2e-4$ H2 for 2.5 min (total: 5 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#6 - 66 degrees
 p_prep = 2.1e-8
 Done 1:30 PM
- H2 Dose#7:** $2e-4$ H2 for 3 min (total: 8 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#7 - 66 degrees
 p_prep = 2.5e-8
 Done at 1:51 PM Turned off source 3 min later 1:54 PM
 Regime of autocatalytic reduction of PdO
- H2 Dose#8:** $2e-4$ H2 for 1 min (total: 9 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#8 - 66 degrees
 p_prep = 3.0e-8
 Done at 2:17 PM
 Regime of autocatalytic reduction of PdO

Adjusted H2 leak valve, pressure now drops to $\sim 2-3e-8$ after 3 min wait after H2 dose and $1.0e-8$ after scan

- H2 Dose#9:** $2e-4$ H2 for 1 min (total: 10 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#9 - 66 degrees

p_prep = 2.0e-8

Done at 2:35 PM

Regime of autocatalytic reduction of PdO

- **H2 Dose#10:** 2e-4 H2 for 5 min (total: 15 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#10 - 66 degrees
p_prep = 1.9e-8
Done at 3:01 PM
- **H2 Dose#11:** 2e-4 H2 for 5 min (total: 20 min) at 298 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#11 - 66 degrees
p_prep = 2.1e-8
Done at 3:24 PM
- **H2 Dose#12:** 2e-4 H2 for 10 min (total: 30 min) at 298 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#12 - 66 degrees
p_prep = 2.2e-8
Done at 3:55 PM
- **H2 Dose#13:** 2e-4 H2 for 10 min (total: 40 min) at 298 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#13 - 66 degrees
p_prep = 2.5e-8
Done at 4:23 PM
Spectra overlapping with #12, but last dose for consistency
- **H2 Dose#14:** 2e-4 H2 for 10 min (total: 50 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#14 - 66 degrees
p_prep = .e-8
Done at PM
Overlapping --> DONE!
Don't forget VB spectrum!!

Background

0.01*5*60=3.0 L per pump-down

13*3 = 39 L per whole day

E-4 H2 Dose equivalent

100 L per second --> 39/100 = 0.39 sec for a E-4 Dose equates accumulated estimates background dose throughout the day in idle period

Legend:

H2 Dose#14 (day before)

After oxidation

H2 Dose#1

H2 Dose#2

H2 Dose#3

H2 Dose#4

H2 Dose#5

H2 Dose#6

H2 Dose#7

H2 Dose#8

H2 Dose#9



H2 Dose#10



H2 Dose#11

H2 Dose#12

H2 Dose#13

H2 Dose#14

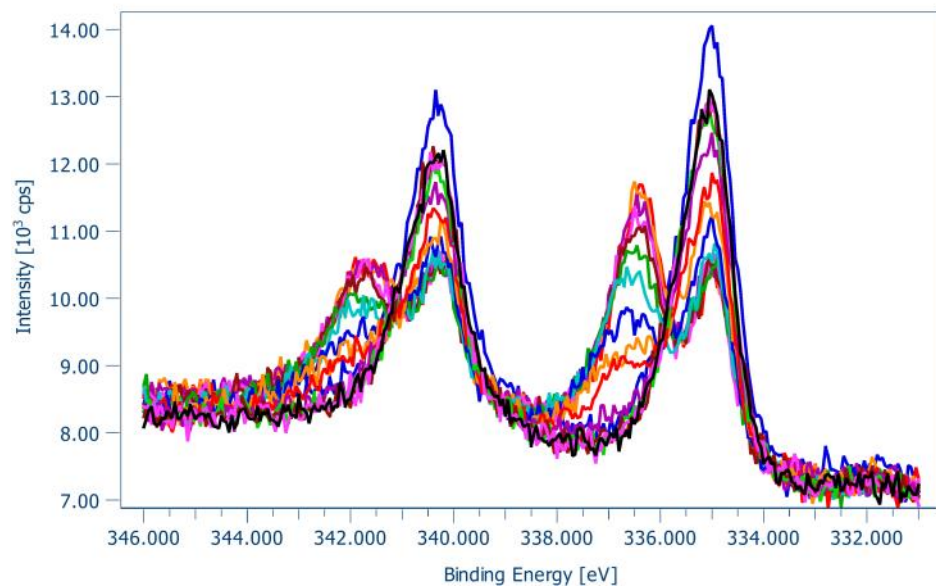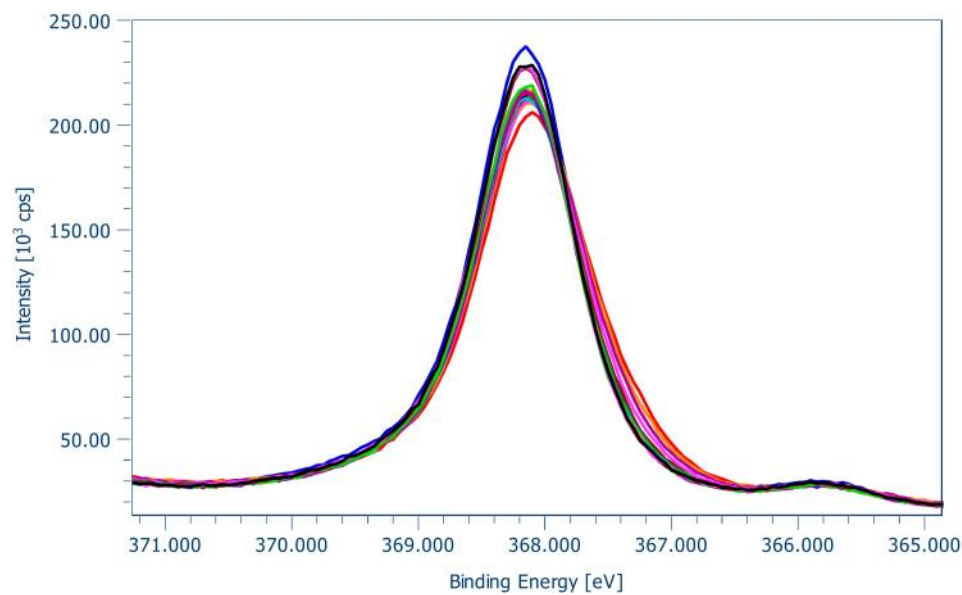

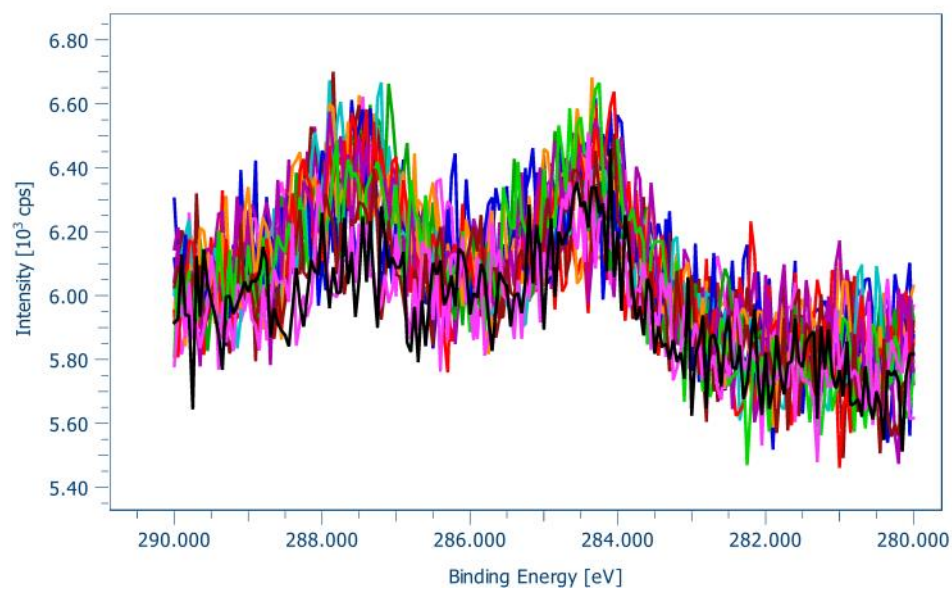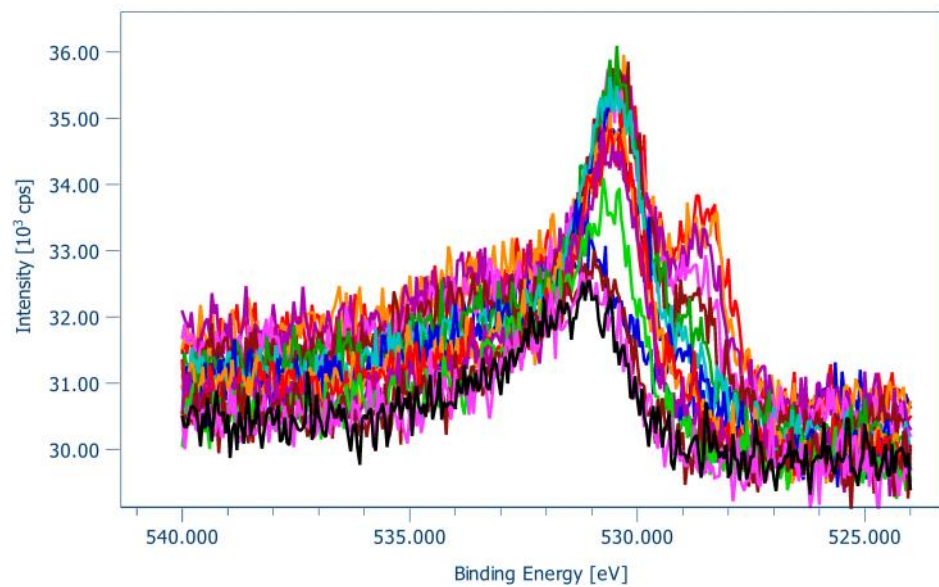

Turned on O2 at 3:50 PM
Turn up O3 at 4:14 PM
Ozone into HP cell 5:18 PM ~30mTorr

200209 - PdAg + 550 K O₂ + H₂ [4th cycle]

Sunday, February 9, 2020 11:59 AM

p_{prep} = 3.6e-9 (background pressure after pumping down over night)

Turn down O₃ and off (12:00 PM)

p_{O₂} = 2025 psi at 12 PM (lost ~375 psi since ~4 PM yday)

Measure

XPS: Ag after 3rd reduction - 66 degrees

Need to leave 1 h 30 mins in UHV (at least till 1:30 PM) between turning down O₃ and oxidizing


Connect cooling

2:02 PM

Turn up O₂ at 2:30 PM

- **Oxidized in 1.92(6) Torr O₂ at 550 K for 5 min**
 - Valve cracked
 - O₂ flowing from O₃ conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24, 25 W to ramp (in 20 sec intervals); let climb to 550 K at 25 W (started timer when it hit 550 K); 21.5 W laser power to stabilize at 550 K
 - Cooling was much stronger this time
 - maintained T within 550-555 K over duration of oxidation
 - Cooling in O₂ till T = 305 K (~1 min), then evacuating HP cell till p_{baratron} = 0.1 mTorr 11:31 AM

Temperature profile recorded with QMS: 200209-file1 (T at beginning 285 K)

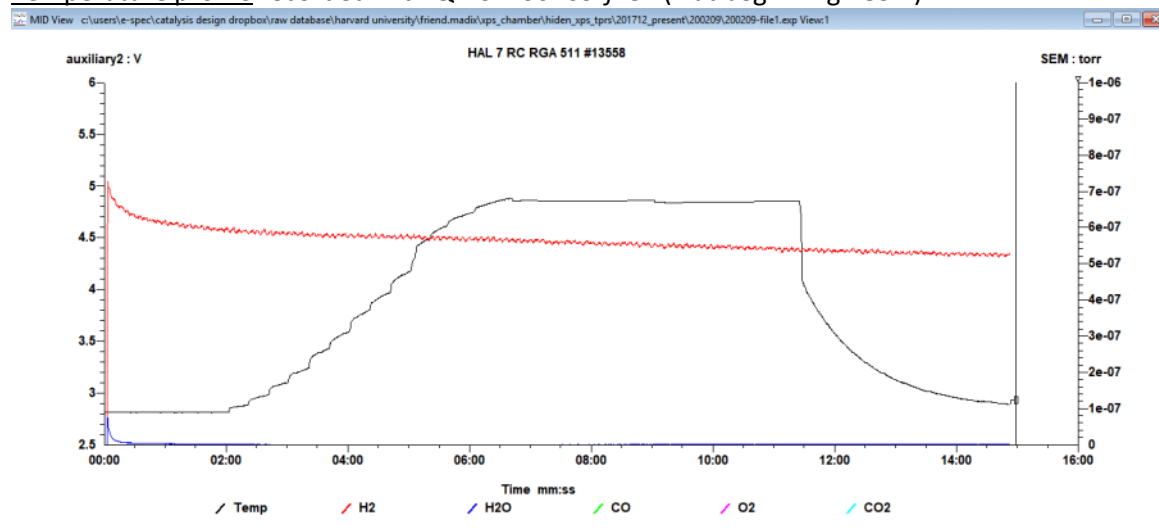

Switch back to air cooling for H₂ dose; clear water from cryostat

Shut valve to HP cell at 15.0"

Turn on X-ray source at 15.0"

O₂ pre-stage pressure was 1975 psi at 3:10 PM

Closed O2 valves



Measure (3:09 PM)


XPS: after oxidation - 66 degrees

H2 Dose#14

Ag after 3rd reduction

After oxidation

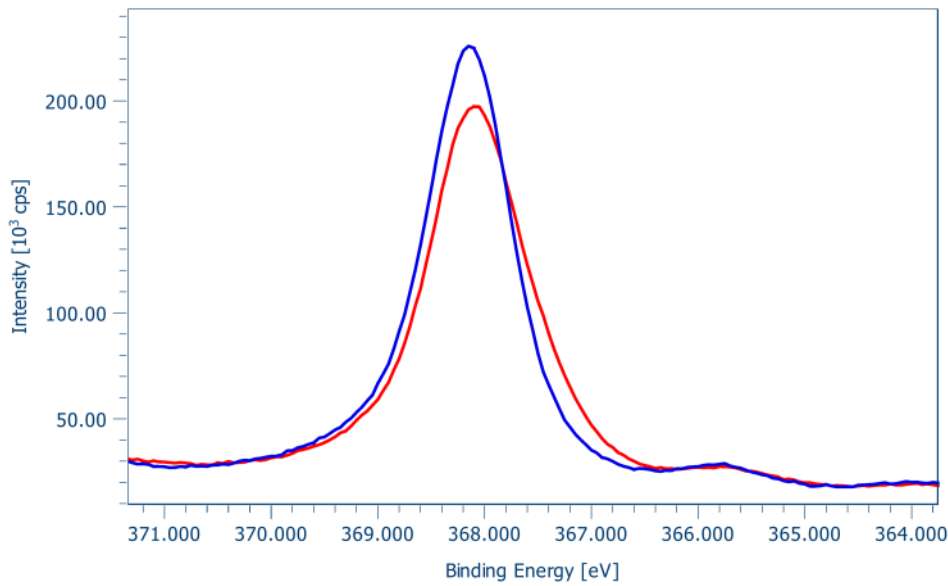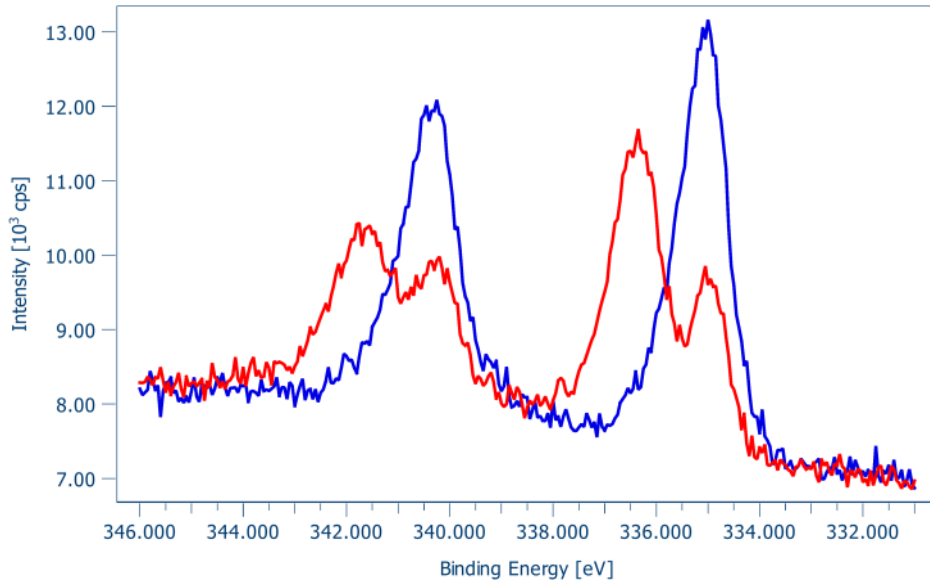

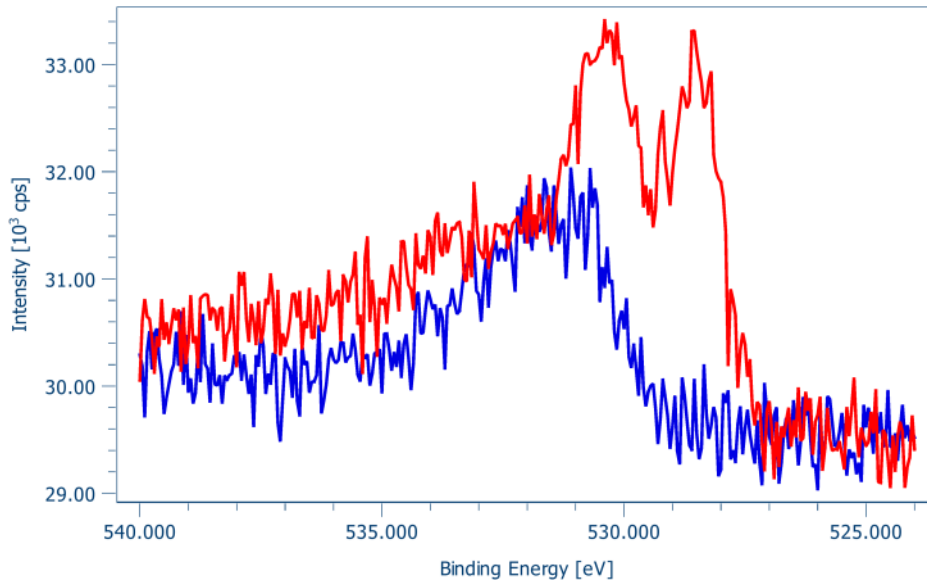

SAME DOSES AS 1ST CYCLE

- H2 Dose#1:** $1e-7$ H2 for 1 min at 295 K (background, right leak valve)
 (4.9e-8 in prep; $p_{\text{bkg}} = 3.6e-9$ $p_{\text{delta}} = 4.5e-8$; had analysis ion gauge off)
XPS: H2 Dose#1 - 66 degrees
 $p_{\text{prep}} = 4.6e-9$
 Done 3:33 PM
- H2 Dose#2:** $1e-5$ H2 for 1 min at 296 K (background, right leak valve)
 (5.0e-6 in prep; had analysis ion gauge off)
XPS: H2 Dose#2 - 66 degrees
 $p_{\text{prep}} = 7.5e-9$
 Done 3:46 PM
- H2 Dose#3:** $1e-4$ H2 for 1 min at 296 K (background, right leak valve)
 (5.0e-5 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure
XPS: H2 Dose#3 - 66 degrees
 $p_{\text{prep}} = 2.0e-8$
 Done 4:02 PM
- H2 Dose#4:** $2e-4$ H2 for 1 min (total: 1 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure
XPS: H2 Dose#4 - 66 degrees
 $p_{\text{prep}} = 2.2e-8$
 Done 4:19 PM
- H2 Dose#5:** $2e-4$ H2 for 1 min (total: 2 min) at 297 K (background, right leak valve)
 (1.0e-4 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure
XPS: H2 Dose#5 - 66 degrees
 $p_{\text{prep}} = 2.0e-8$
 Done 4:37 PM

- **H2 Dose#6:** 2e-4 H2 for 3 min (total: 5 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#6 - 66 degrees
 p_prep = 2.1e-8
 Done 4:58 PM
- **H2 Dose#7:** 2e-4 H2 for 3 min (total: 8 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#7 - 66 degrees
 p_prep = 2.2e-8
 Done at PM
- **H2 Dose#8:** 2e-4 H2 for 1 min (total: 9 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#8 - 66 degrees
 p_prep = 2.5e-8
 Done at 5:39 PM
- **H2 Dose#9:** 2e-4 H2 for 1 min (total: 10 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#9 - 66 degrees
 p_prep = 1.8e-8
 Done at 5:59 PM
 Regime of autocatalytic reduction of PdO?
- **H2 Dose#10:** 2e-4 H2 for 5 min (total: 15 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#10 - 66 degrees
 p_prep = 2.1e-8
 Done at 6:24 PM

Probably start O2 flow again
- **H2 Dose#11:** 2e-4 H2 for 5 min (total: 20 min) at 298 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#11 - 66 degrees
 p_prep = 1.6e-8
 Done at 6:49 PM
- **H2 Dose#12:** 2e-4 H2 for 10 min (total: 30 min) at 298 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
 Wait 3 min to recover pressure
 XPS: H2 Dose#12 - 66 degrees
 p_prep = 2.1e-8
 Done at 7:17 PM
- **H2 Dose#13:** 2e-4 H2 for 10 min (total: 40 min) at 298 K (background, right leak valve)

(1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#13 - 66 degrees

p_prep = 2.2e-8

Done at 4:23 PM

Spectra overlapping with #12, but last dose for consistency

- **H2 Dose#14:** 2e-4 H2 for 10 min (total: 50 min) at 297 K (background, right leak valve)
(1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#14 - 66 degrees

p_prep = 2.0e-8

Done at PM

Overlapping --> DONE!

Don't forget VB spectrum!!

Legend:

~~H2 Dose#14 (day before)~~

Ag after 3rd reduction

After oxidation

H2 Dose#1

H2 Dose#2

H2 Dose#3

H2 Dose#4

H2 Dose#5

H2 Dose#6

H2 Dose#7

H2 Dose#8

H2 Dose#9


H2 Dose#10




H2 Dose#11

H2 Dose#12

H2 Dose#13

H2 Dose#14

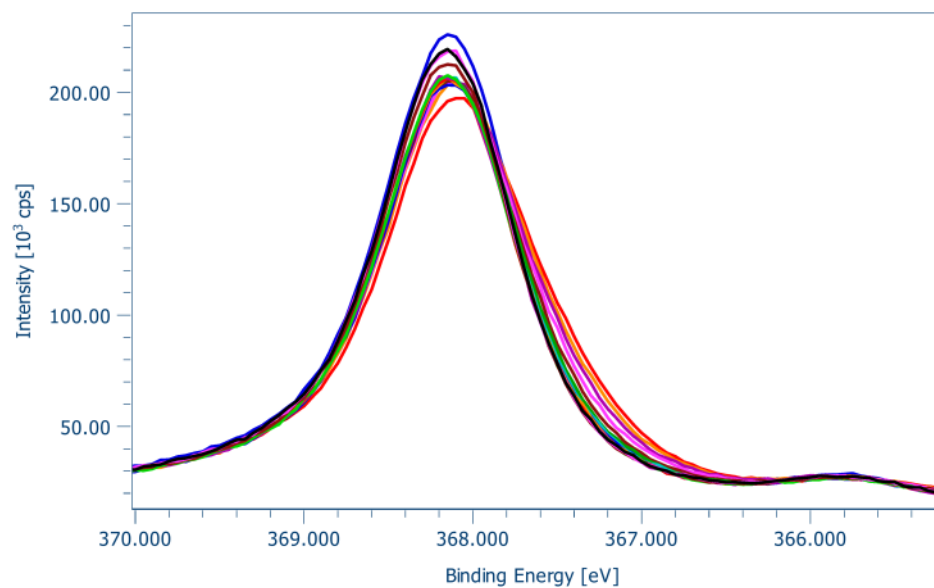

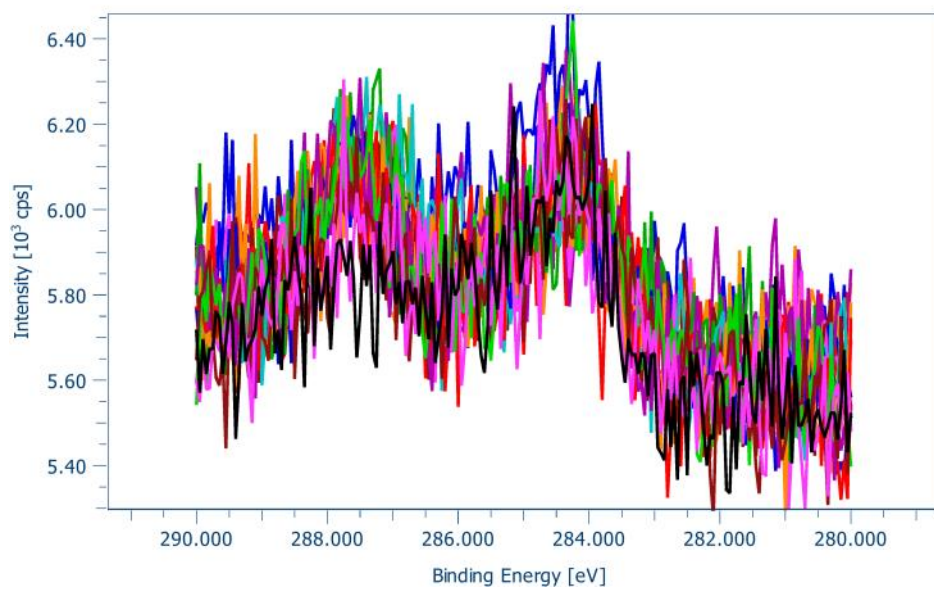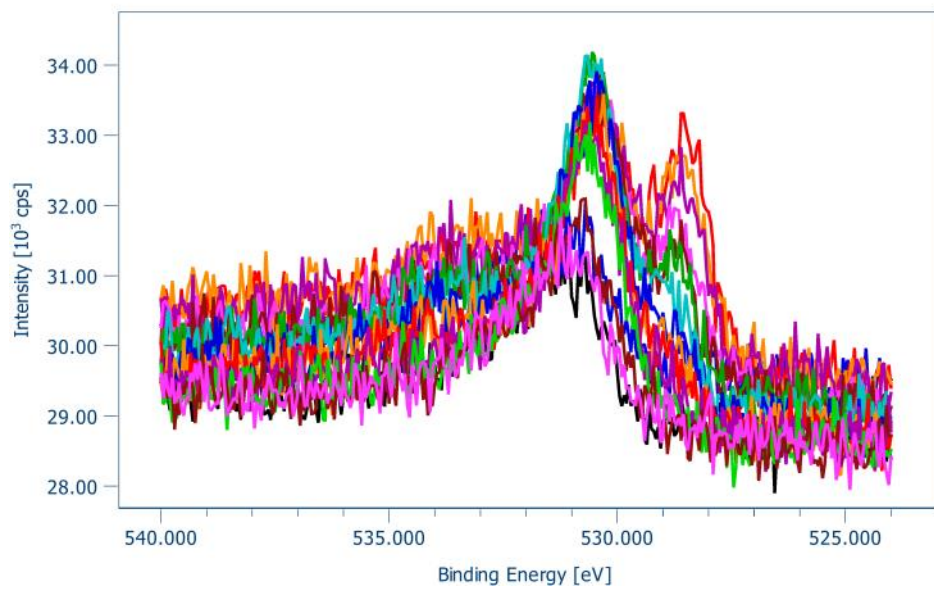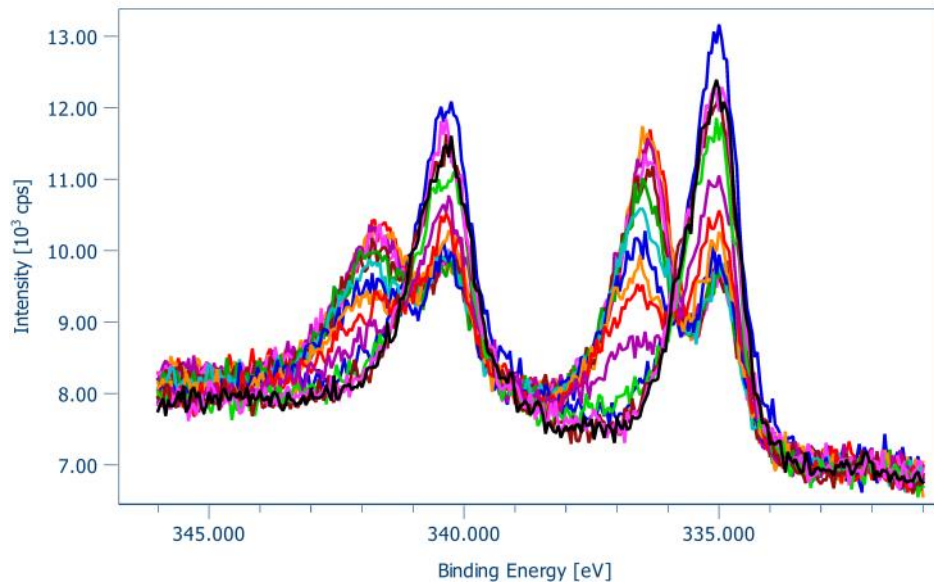

Turned on O2 at 5:50 PM

Turn up O3 at 6:25 PM

Ozone into HP cell 8:25 PM ~30mTorr

200210 - PdAg + 550 K O₂ + H₂ [5th cycle]

Monday, February 10, 2020 9:19 AM

p_{prep} = 2.4e-9 (background pressure after pumping down over night)

Turn down O₃ and off (9:20 AM)

p_{O₂} = 1750 psi at 9:20 AM (lost ~225 psi since ~9 PM yday)

Measure

XPS: Ag after 4th reduction - 66 degrees

Need to leave 1 h 30 mins in UHV (at least till 10:50 AM) between turning down O₃ and oxidizing


Connect cooling

10:30 AM

Turn up O₂ at 11:24 AM

- **Oxidized in 1.93(1) Torr O₂ at 550 K for 5 min**
 - Valve cracked
 - O₂ flowing from O₃ conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24, 25 W to ramp (in 20 sec intervals); let climb to 550 K at 25 W (started timer when it hit 550 K); 23.5 W laser power to stabilize at 550 K
 - maintained T within 550-555 K over duration of oxidation
 - Cooling in O₂ till T = 305 K (~1 min) , then evacuating HP cell till p_{baratron} = 0.1 mTorr 11:50 AM

Temperature profile recorded with QMS: 200210-file1 (T at beginning 286 K)

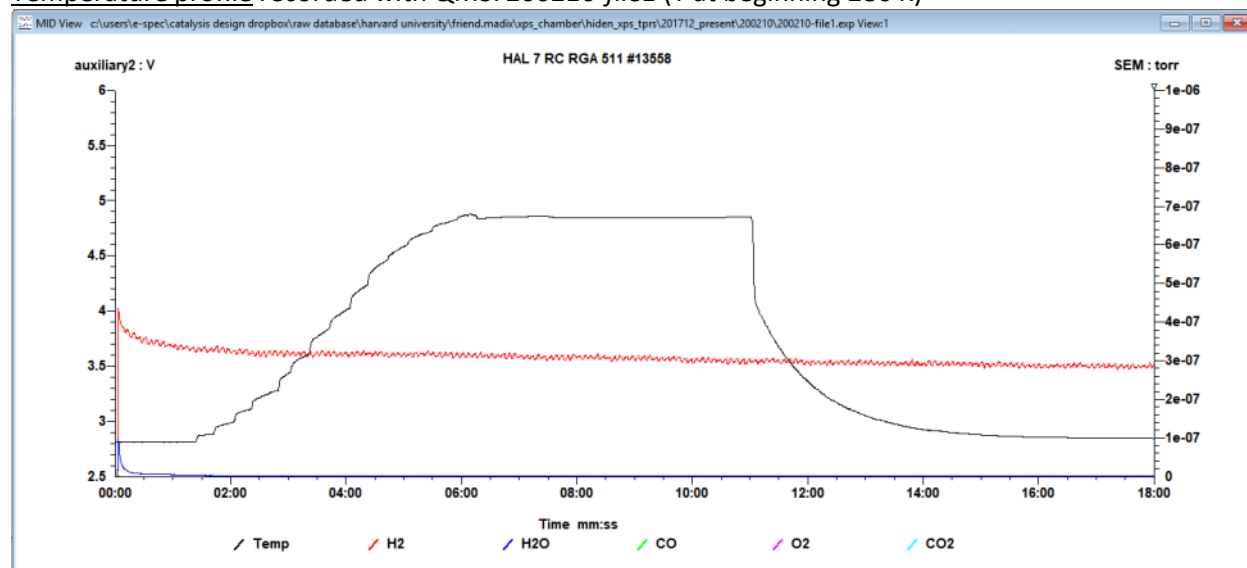

Switch back to air cooling for H₂ dose; clear water from cryostat

Shut valve to HP cell at 15.0"

Turn on X-ray source at 15.0"

O₂ pre-stage pressure was 1700 psi at 12:04 PM

Closed O2 valves



Measure (12:01 PM)



XPS: after oxidation - 66 degrees

H2 Dose#14

Ag after 2nd reduction

After oxidation

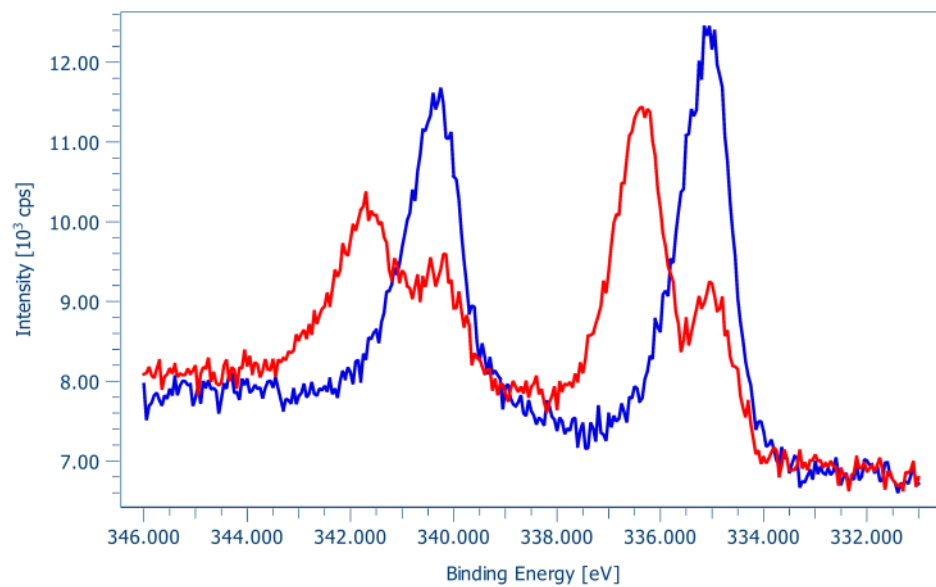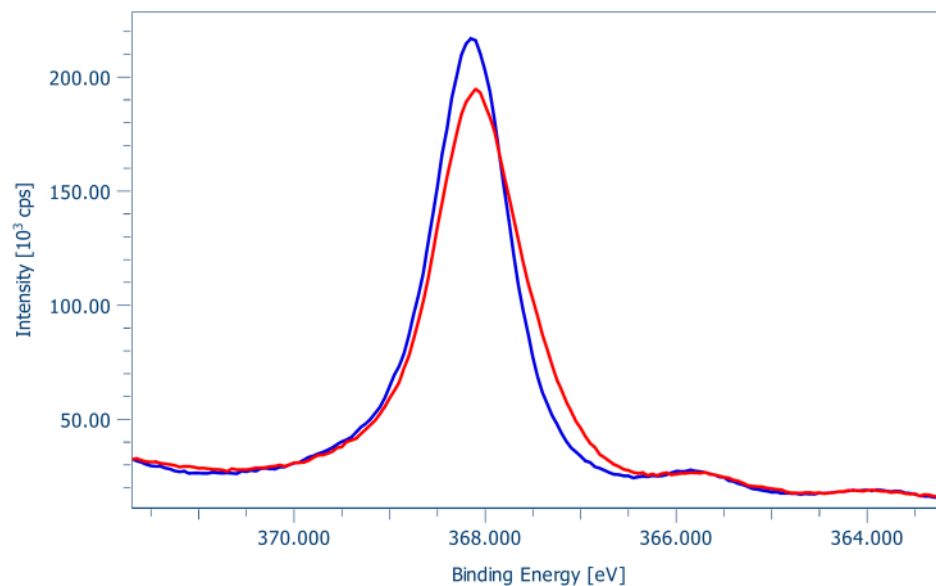

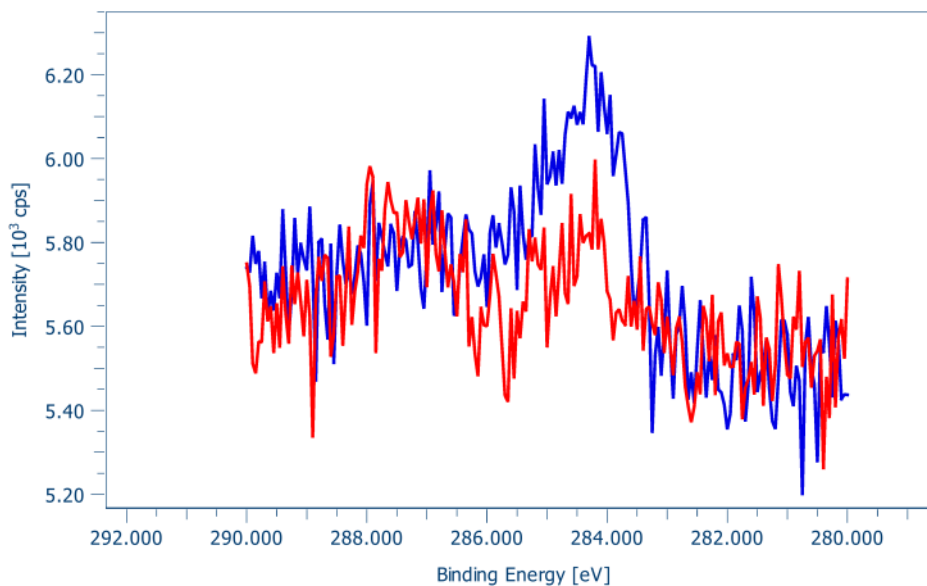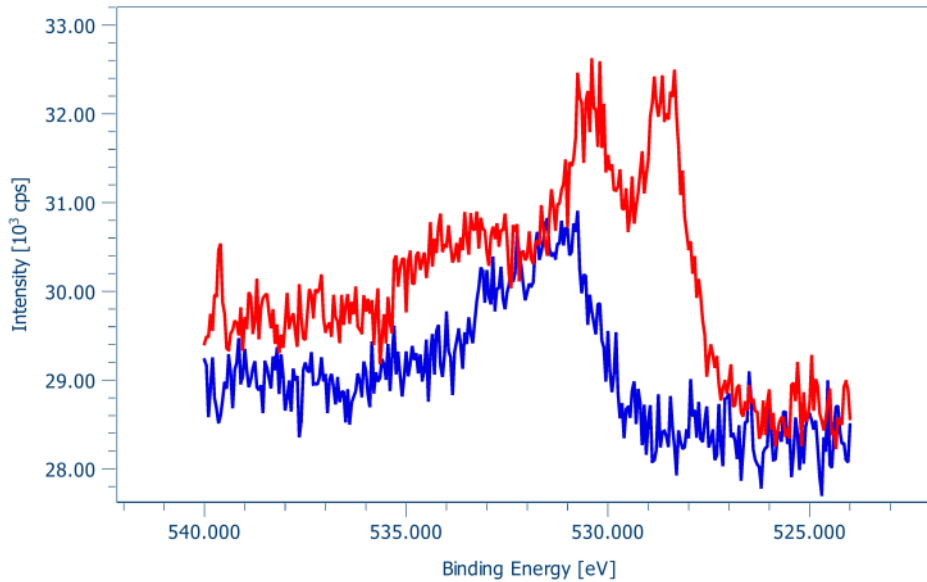

SAME DOSES AS 1ST CYCLE

- H2 Dose#1:** $1e-7$ H2 for 1 min at 295 K (background, right leak valve)
 (4.8e-8 in prep; $p_{\text{bkg}} = 2.7e-9$ $p_{\text{delta}} = 4.5e-8$; had analysis ion gauge off)
XPS: H2 Dose#1 - 66 degrees
 $p_{\text{prep}} = 3.8e-9$
 Done 12:23 PM
- H2 Dose#2:** $1e-5$ H2 for 1 min at 296 K (background, right leak valve)
 (5.0e-6 in prep; had analysis ion gauge off)
XPS: H2 Dose#2 - 66 degrees
 $p_{\text{prep}} = 6.4e-9$
 Done 12:38 PM
- H2 Dose#3:** $1e-4$ H2 for 1 min at 296 K (background, right leak valve)
 (5.0e-5 in prep; had analysis ion gauge off)
 Wait 2.5 min to recover pressure

XPS: H2 Dose#3 - 66 degrees

p_prep = 1.9e-8

Done 12:54 PM

- **H2 Dose#4:** 2e-4 H2 for 1 min (total: 1 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 2.5 min to recover pressure

XPS: H2 Dose#4 - 66 degrees

p_prep = 1.9e-8

Done 1:11 PM

- **H2 Dose#5:** 2e-4 H2 for 1 min (total: 2 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 2.5 min to recover pressure

XPS: H2 Dose#5 - 66 degrees

p_prep = 2.1e-8

Done 1:28 PM

- **H2 Dose#6:** 2e-4 H2 for 3 min (total: 5 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#6 - 66 degrees

p_prep = 2.2e-8

Done 1:47 PM

- **H2 Dose#7:** 2e-4 H2 for 3 min (total: 8 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#7 - 66 degrees

p_prep = 1.8e-8

Done at 2:06 PM

- **H2 Dose#8:** 2e-4 H2 for 1 min (total: 9 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#8 - 66 degrees

p_prep = 2.0e-8

Done at 2:23 PM

- **H2 Dose#9:** 2e-4 H2 for 1 min (total: 10 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#9 - 66 degrees

p_prep = 1.9e-8

Done at 2:39 PM

~Regime of autocatalytic reduction of PdO

- **H2 Dose#10:** 2e-4 H2 for 5 min (total: 15 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)

Wait 3 min to recover pressure

XPS: H2 Dose#10 - 66 degrees

p_prep = 1.9 e-8

Done at 3:01 PM

~Regime of autocatalytic reduction of PdO

- **H2 Dose#11:** 2e-4 H2 for 5 min (total: 20 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#11 - 66 degrees
p_prep = 2.1e-8
Done at 3:23 PM
- **H2 Dose#12:** 2e-4 H2 for 10 min (total: 30 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#12 - 66 degrees
p_prep = 2.1e-8
Done at 3:51 PM
- **H2 Dose#13:** 2e-4 H2 for 10 min (total: 40 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#13 - 66 degrees
p_prep = 2.2e-8
Done at 4:19 PM
- **H2 Dose#14:** 2e-4 H2 for 10 min (total: 50 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#14 - 66 degrees
p_prep = 1.7e-8
Done at 4:51 PM
Background lower for some reason (what happened??), do one more dose and measure to check if spectra are overlapping
Don't forget VB spectrum!!
- **H2 Dose#15:** 2e-4 H2 for 10 min (total: 60 min) at 297 K (background, right leak valve) (1.0e-4 in prep; had analysis ion gauge off)
Wait 3 min to recover pressure
XPS: H2 Dose#15 - 66 degrees
p_prep = 2.3e-8
Done at PM
Basically overlapping spectra compared to H2 Dose#13 but still small baseline offset (more visible at higher BE e.g. for O1s)
Something with H2 Dose#14 wasn't right and the baseline was offset. Should probably use Dose#15 and replace with #14 (total dosing time of 50 min)
Don't forget VB spectrum!!

Legend:

~~H2 Dose#14 (day before)~~

Ag after 4th reduction



After oxidation



H2 Dose#1

H2 Dose#2

H2 Dose#3

- H2 Dose#4
- H2 Dose#5
- H2 Dose#6
- H2 Dose#7
- H2 Dose#8
- H2 Dose#9
- H2 Dose#10
- H2 Dose#11
- H2 Dose#12
- H2 Dose#13
- H2 Dose#14
- H2 Dose#15

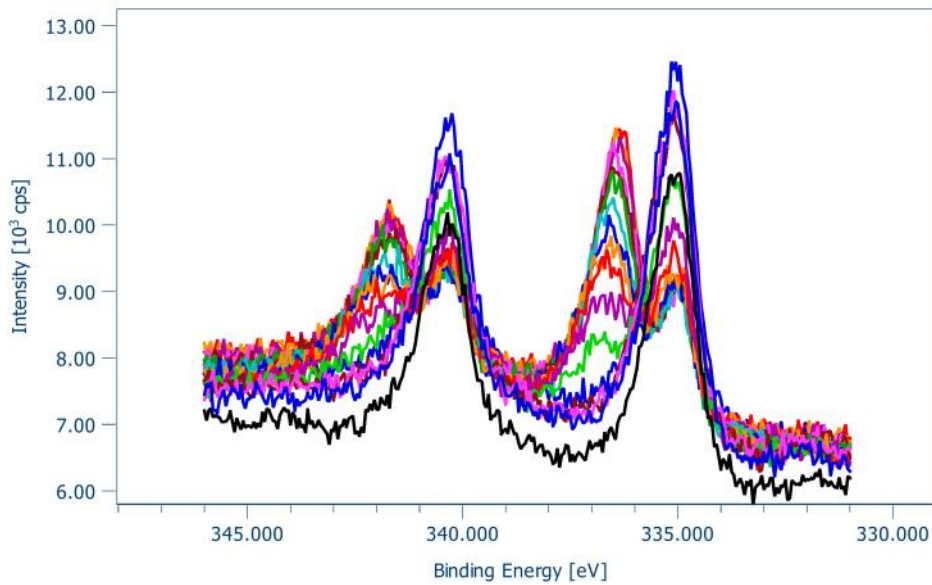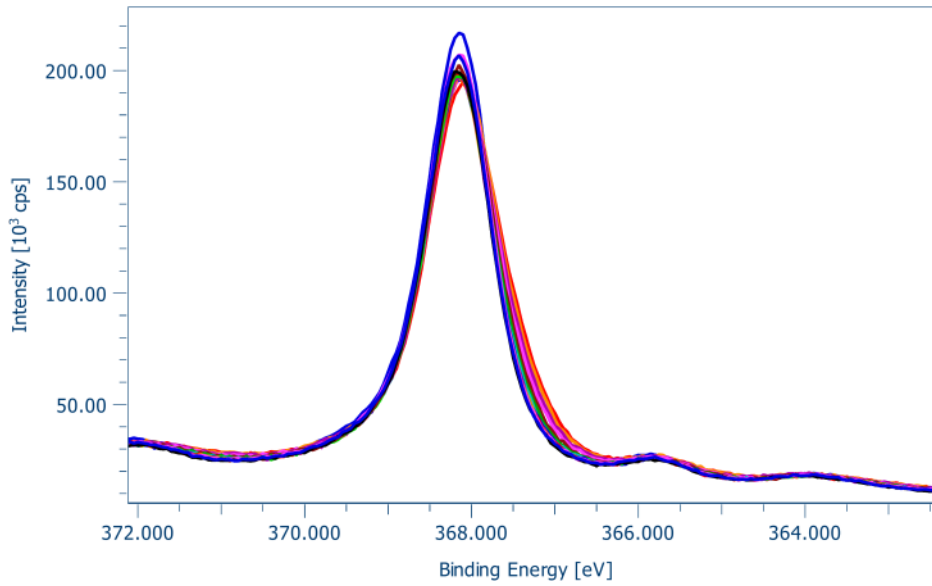

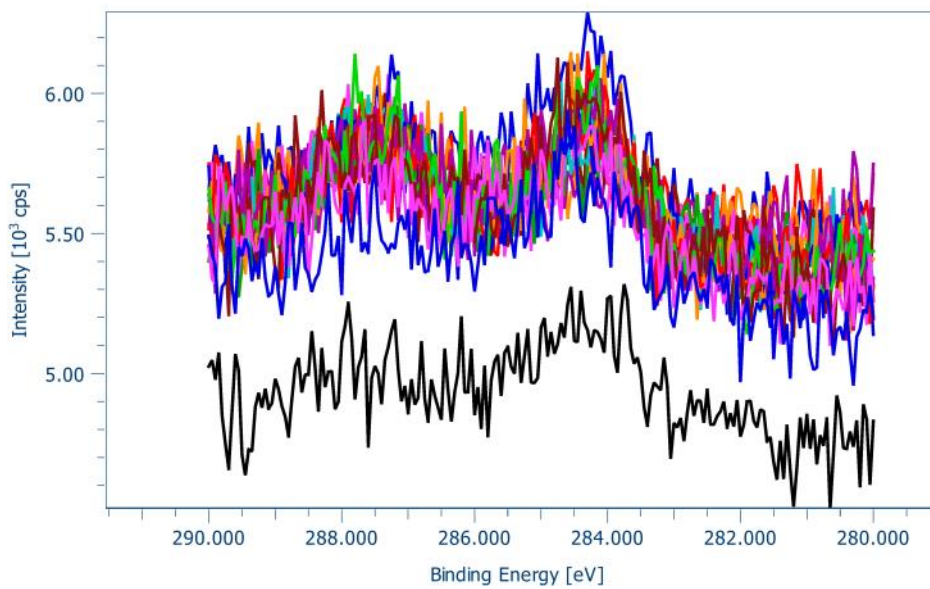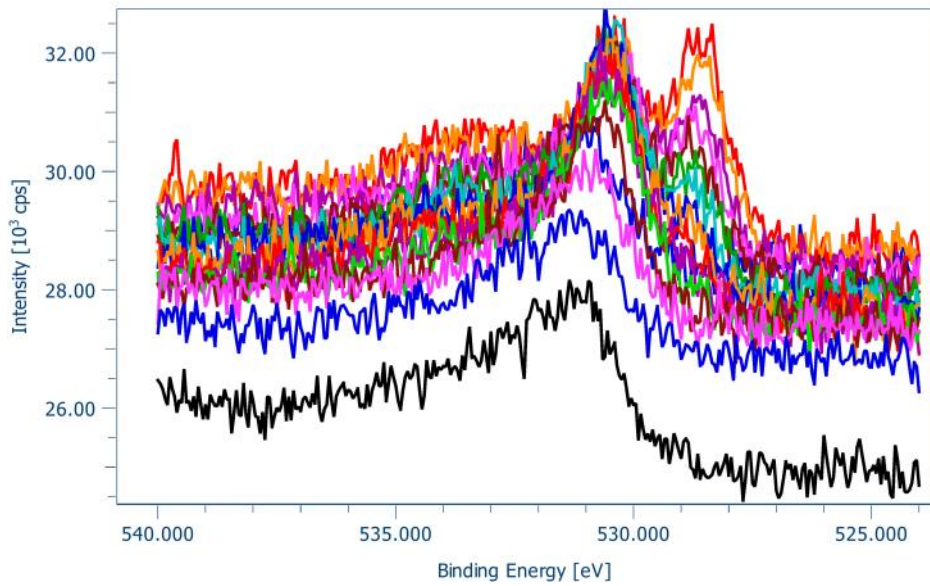

Turned on O2 at 4:22 PM
Turn up O3 at 4:55 PM
Ozone into HP cell 6:15 PM ~30mTorr

200218 - Ag + 550 K O₂ [control]

Tuesday, February 18, 2020 9:51 AM


Turn down O₃ and off (9:30 AM)

p_{O₂} = 1400 psi at 10:00 AM

Measure

XPS: Ag after sputter - 66 degrees

Detector must have suffered some loss over the past experiments :/

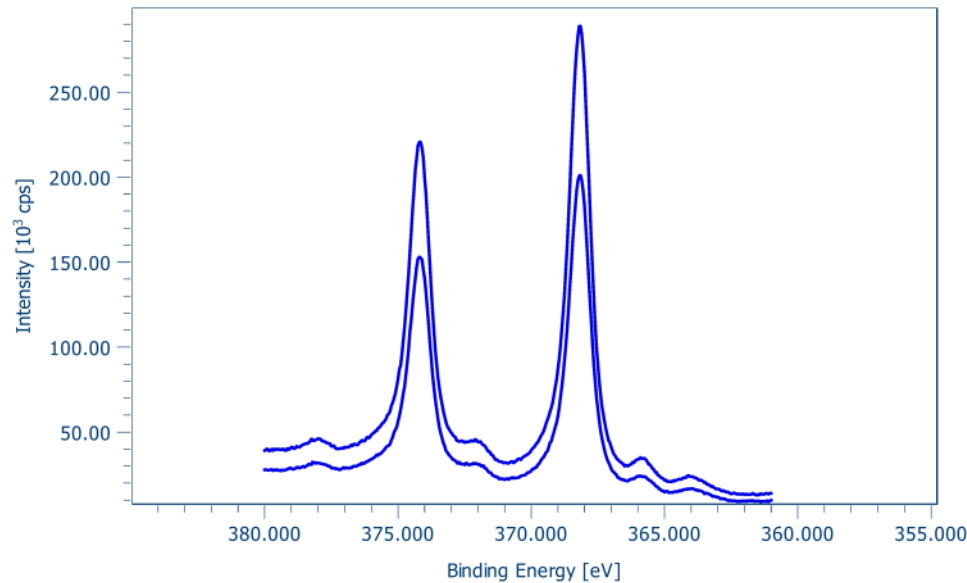

Need to leave 1 h 30 mins in UHV (at least till 11:00 AM) between turning down O₃ and oxidizing


Connect cooling

11:10 AM

Turn up O₂ at 11:20 AM

- **Oxidized in 1.93(0) Torr O₂ at 550 K for 5 min**
 - Valve cracked
 - O₂ flowing from O₃ conditions before
 - turbo 19 W
 - 1, 3, 5, 7, 10, 12, 15, 17, 20, 23, 24, 25, 26, 27, 28 W to ramp (in 20 sec intervals); let climb to 550 K at 28 W (started timer when it hit 550 K); 27 W laser power to stabilize at 550 K
 - Cooling pretty strong this time, therefore higher laser power
 - maintained T within 550-555 K over duration of oxidation
 - Cooling in O₂ till T = 305 K (~1 min), then evacuating HP cell till p_{baratron} = 0.1 mTorr 11: AM

Temperature profile recorded with QMS: 200218-file1 (T at beginning 284 K)

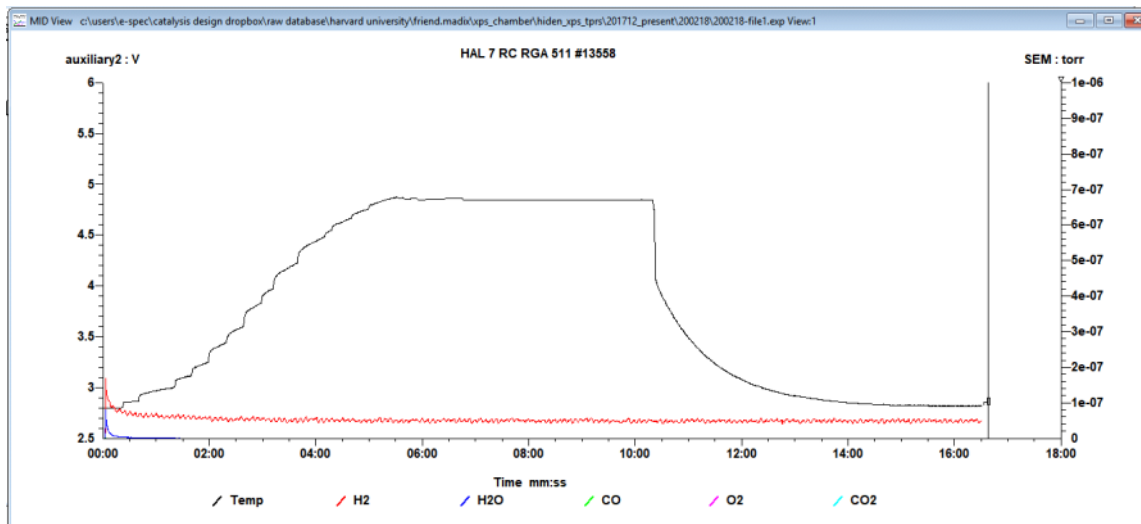

Switch back to air cooling for H2 dose; clear water from cryostat
 Shut valve to HP cell at 15.0"
 Turn on X-ray source at 15.0"


O2 pre-stage pressure was 1350 psi at 12:00 PM
 Closed O2 valves


Measure (12:00 PM)

XPS: after oxidation - 66 degrees

Ag after sputter

After oxidation

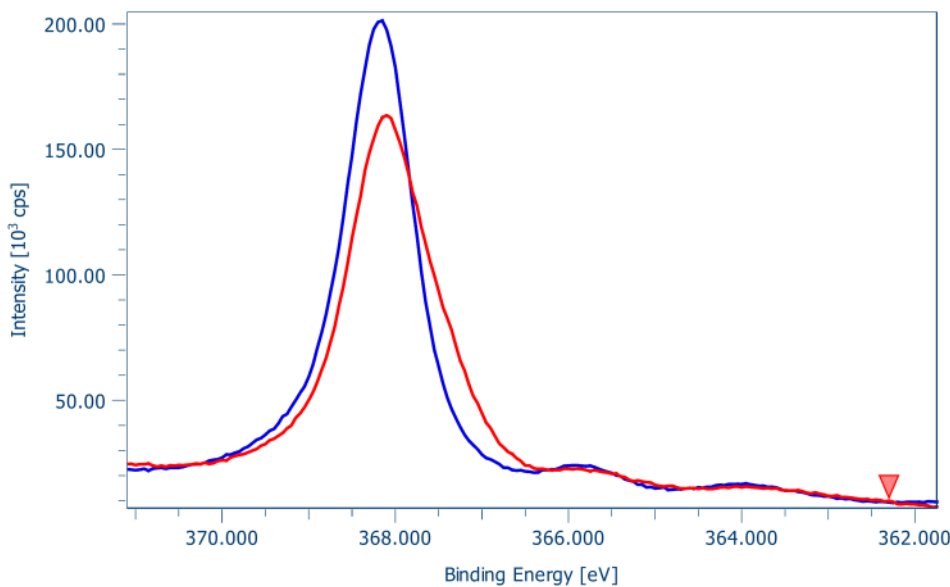

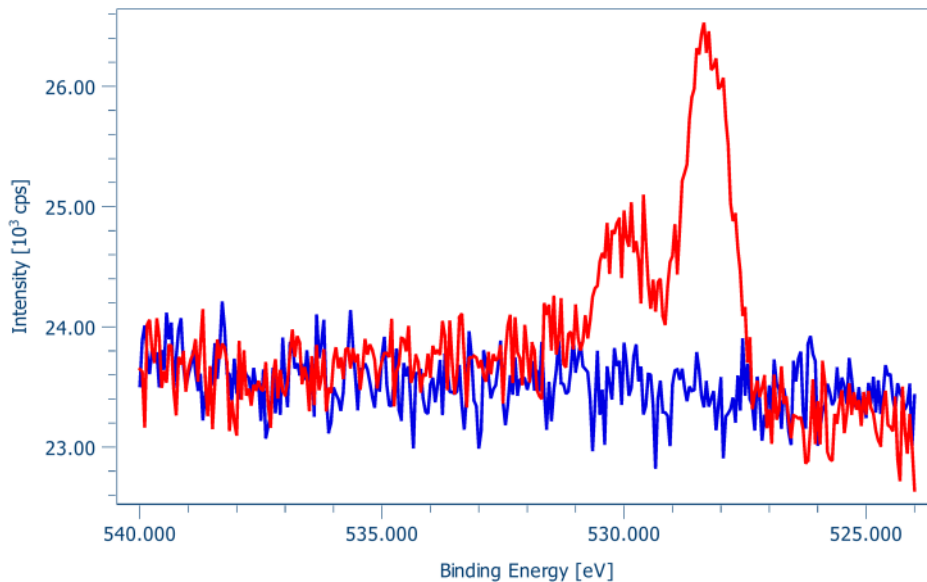

SAME DOSES AS 1ST CYCLE but don't dose H2

- **H2 Dose#1:** 1 min at 295 K
(had analysis ion gauge off)
XPS: H2 Dose#1 - 66 degrees
p_prep = 1.6e-9
Done 12:22 PM

- **H2 Dose#2:** 1 min at 296 K
(had analysis ion gauge off)
XPS: H2 Dose#2 - 66 degrees
p_prep = 1.5e-9
Done 12:34 PM

- **H2 Dose#3:** 1 min at 296 K
(had analysis ion gauge off)
XPS: H2 Dose#3 - 66 degrees
p_prep = 1.5e-9
Done 12:45 PM

- **H2 Dose#4:** 1 min (total: 1 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#4 - 66 degrees
p_prep = 1.5e-9
Done 12:57 PM

- **H2 Dose#5:** 1 min (total: 2 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#5 - 66 degrees
p_prep = 1.5e-9
Done 1:08 PM

- **H2 Dose#6:** 3 min (total: 5 min) at 297 K
(had analysis ion gauge off)

XPS: H2 Dose#6 - 66 degrees

p_prep = 1.4e-9

Done 1:25 PM

- **H2 Dose#7:** 3 min (total: 8 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#7 - 66 degrees
p_prep = 1.3e-9
Done at 1:40 PM
- **H2 Dose#8:** 1 min (total: 9 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#8 - 66 degrees
p_prep = 1.4e-9
Done at 1:53 PM
- **H2 Dose#9:** 1 min (total: 10 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#9 - 66 degrees
p_prep = 1.3e-9
Done at 2:06 PM
- **H2 Dose#10:** 5 min (total: 15 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#10 - 66 degrees
p_prep = 1.3e-9
Done at 2:22 PM
- **H2 Dose#11:** 5 min (total: 20 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#11 - 66 degrees
p_prep = 1.3e-9
Done at 2:40 PM
- **H2 Dose#12:** 10 min (total: 30 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#12 - 66 degrees
p_prep = 1.3e-9
Done at 3:05 PM
- **H2 Dose#13:** 10 min (total: 40 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#13 - 66 degrees
p_prep = 1.3e-9
Done at 3:29 PM
- **H2 Dose#14:** 10 min (total: 50 min) at 297 K
(had analysis ion gauge off)
XPS: H2 Dose#14 - 66 degrees
p_prep = 1.3e-9
Done at 3:55 PM
- **H2 Dose#15:** 10 min (total: 50 min) at 297 K

(had analysis ion gauge off)
XPS: H2 Dose#15 - 66 degrees
p_prep = 1.3e-9
Done at 4:22 PM

Legend:

Ag after sputter

After oxidation

H2 Dose#1

H2 Dose#2

H2 Dose#3

H2 Dose#4

H2 Dose#5

H2 Dose#6

H2 Dose#7

H2 Dose#8

H2 Dose#9

H2 Dose#10


H2 Dose#11



H2 Dose#12

H2 Dose#13

H2 Dose#14

H2 Dose#15

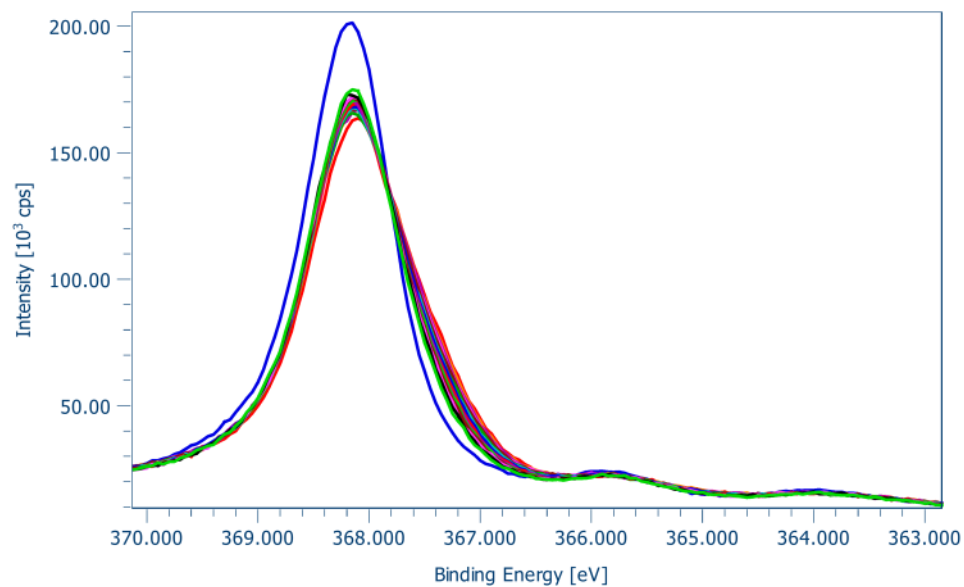

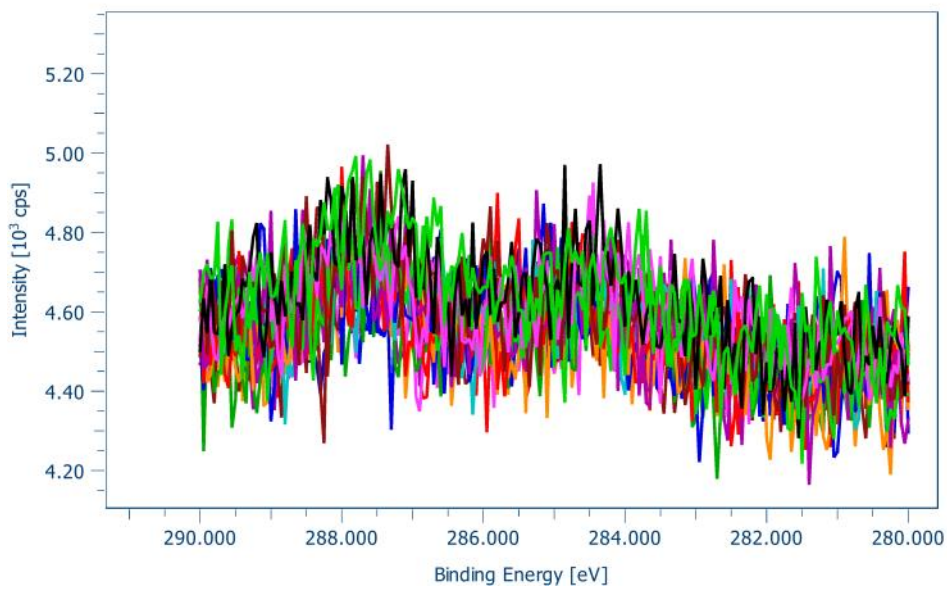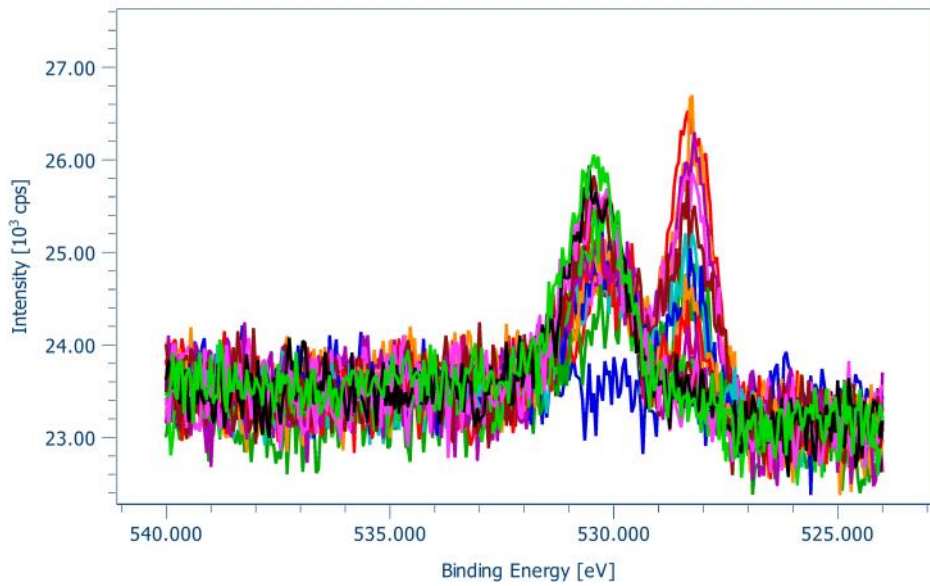

- **TPD:** 200218-file2
 Some CO₂ around 400 K --> from Ag carbonate
 Rest is contribution from Cu / sample holder it seems at > 600 K

